


Treatment Concepts for  
**EXTRACTION  
SOCKETS**



**NEW**

New clinical cases  
Latest scientific literature  
QR codes linked to clinical  
videos, 3D-animation...

# ALL OUR SCIENCE IN YOUR HANDS

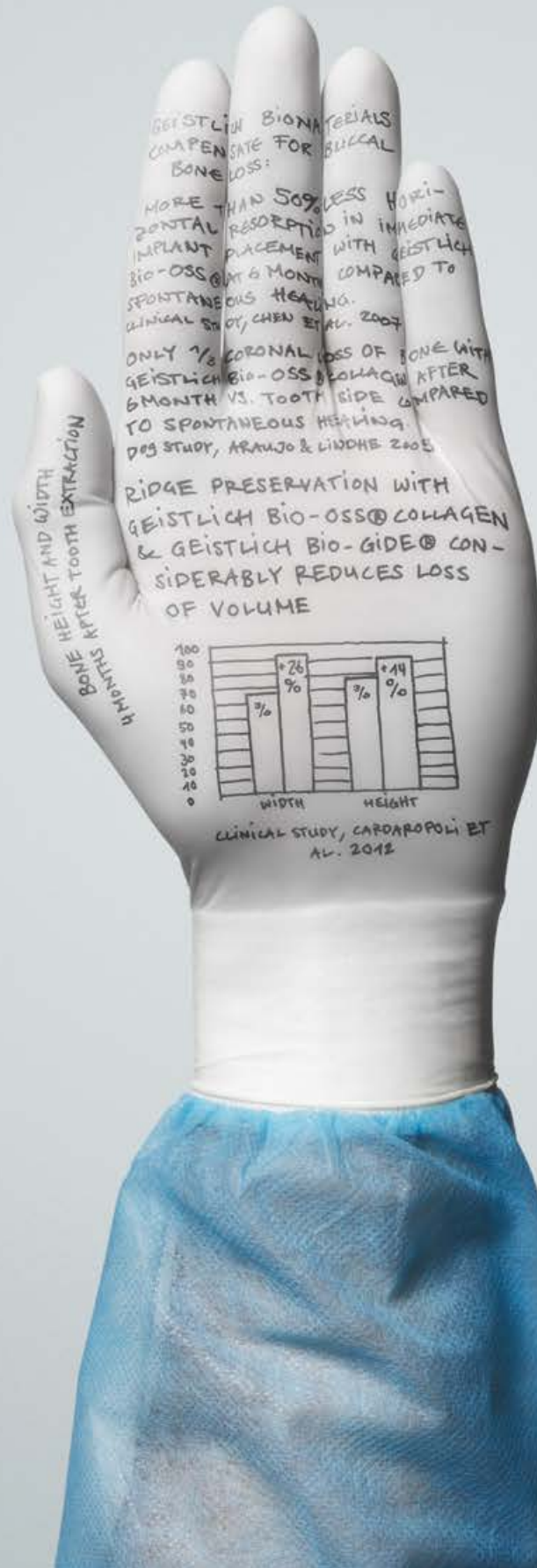


# TABLE OF CONTENTS

3	<b>Clinical evidence</b>	
7	<b>Immediate implant placement</b>	
	Dr. Michael Back & Dr. Oliver Blume (Munich, Germany)	
	Dr. Franck Bonnet (Le Cannet, France)	
11	<b>Early implant placement</b>	
	Prof. Daniel Buser & Prof. Urs Belser (Berne, Switzerland)	
	Dr. Luca Cordaro (Rome, Italy)	
	Dr. Raffaele Cavalcanti (Bari, Italy)	
17	<b>Delayed/late implant placement</b>	
	Dr. Hadi Antoun & Dr. Bouchra Sojod (Paris, France)	
	Dr. Stefan Fickl (Würzburg, Germany)	
21	<b>Extraction socket treatment options</b>	
23	<b>Delayed/late implant placement</b>	
	Dr. Célia Coutinho Alves (Porto, Portugal)	
	Prof. Ronald E. Jung (Zurich, Switzerland)	
	Prof. Carlo Maiorana (Milan, Italy)	
	Prof. Julio Cesar Joly, Prof. Robert Carvalho da Silva &	
	Prof. Paulo Fernando M. de Carvalho (Sao Paulo, Brazil)	
	Dr. Fernán López (Medellin, Colombia)	
	Dr. Ham Byung-Do (Kainos Dental Clinic, Seoul, Korea)	
35	<b>No implant placement</b>	
	Dr. Jeffrey Ganeles (Boca Raton, USA)	
	Dr. Philipp Grohmann (Berikon, Switzerland)	
39	<b>Product Information</b>	

# AFTER TOOTH EXTRACTION:

Spontaneous healing implies

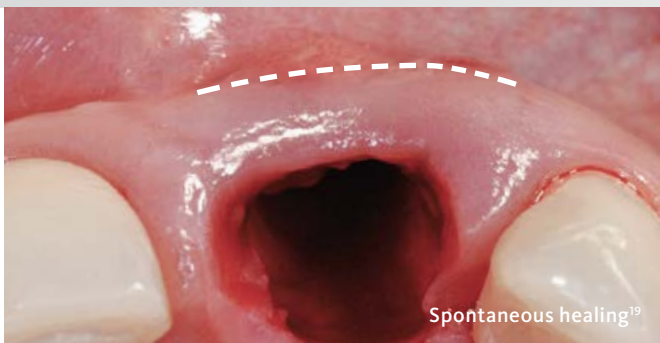

alveolar volume loss<sup>1-5</sup>



## What happens with spontaneous healing?

The healing of extraction sockets and the resorption processes that take place after tooth extraction have been investigated thoroughly in recent years. The most recent clinical studies have shown that:

- › the alveolar volume loss after tooth extraction is severe<sup>1-5</sup>
- › two-thirds of resorption take place within the first three months<sup>1</sup>

## Ridge volume loss after extraction in numbers:

**Horizontal loss:**

**– 49 %<sup>1</sup>**  
(after 12 months)

**– 3.8 mm<sup>4</sup>**  
(after 6 months)

**Vertical loss:**

**from – 1.2 mm<sup>4</sup>**  
(after 6 months)

**to – 1.5 mm<sup>7</sup>**  
(after ca. 6 months)

## Volume loss: clinical implications

**Potentially important clinical implications of spontaneous healing compared to Ridge Preservation:**


- › poorer maintenance of healthy periimplant soft tissues<sup>6</sup>
- › poorer aesthetic outcomes<sup>6</sup>
- › 10 times greater need for hard tissue augmentation at implant placement without previous Ridge Preservation<sup>7</sup>




Systematic review quote:

“We found that alveolar ridge preservation is effective in limiting physiologic ridge reduction as compared with tooth extraction alone.”<sup>8</sup>

Ridge Preservation with Geistlich Biomaterials

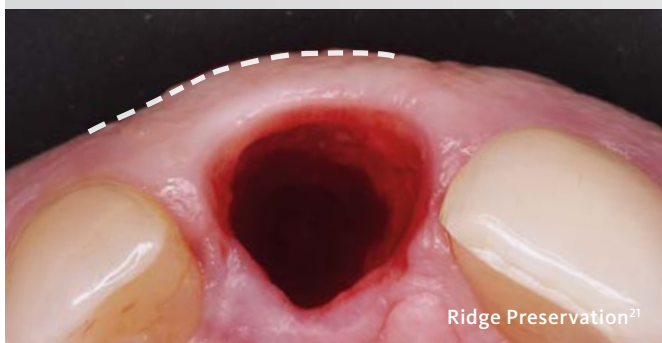

largely maintains the alveolar volume<sup>5,10,11</sup>

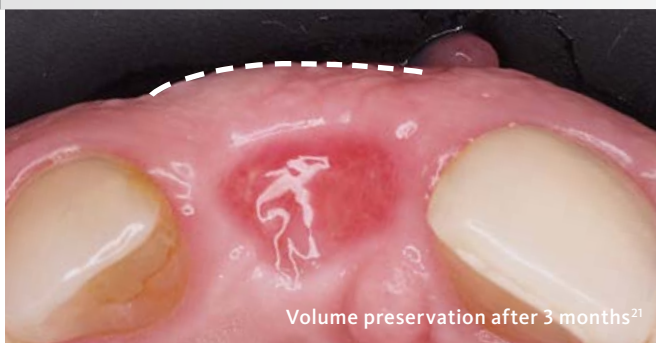

### Ridge Preservation pays off.

While immediate implant placement does not prevent bone resorption<sup>9</sup>, the treating extraction sockets with Geistlich Biomaterials can largely compensate for bone loss and preserve the contour of the alveolar ridge.<sup>5,10,11</sup>

### Volume preservation: clinical evidence

The latest systematic reviews (high level of clinical evidence) agree that Ridge Preservation is effective in limiting alveolar volume loss.<sup>8,12-16</sup>

#### Ridge Preservation with Geistlich Biomaterials can:

- › prevent volume loss and lead to an optimised hard and soft tissue situation irrespective of the chosen time for implantation<sup>17</sup>
- › improve the aesthetic outcome by preserving the alveolar ridge volume and contour, when the objective of treatment is to place a bridge<sup>18</sup>


#### References

- 1 Schropp L, et al. *Int J Periodontics Restorative Dent.* 2003 Aug;23(4):313-23.
- 2 Van der Weijden F, et al. *J Clin Periodontol.* 2009 Dec;36(12):1048-58.
- 3 Sanz M, et al. *Clin Oral Implants Res.* 2010 Jan;21(1):13-21.
- 4 Hämmerle CH, et al. *Clin Oral Implants Res.* 2012 Feb;23 Suppl 5:80-2.
- 5 Jung RE, et al. *J Clin Periodontol.* 2013 Jan;40(1):90-8.
- 6 Vignoletti F, et al. *Clin Oral Implants Res.* 2012 Feb;23 Suppl 5:22-38.
- 7 Weng D, et al. *Eur J Oral Implantol.* 2011;4 Suppl:59-66.
- 8 Avila-Ortiz G, et al. *J Dent Res.* 2014 Oct;93(10):950-8.
- 9 Wang RE & Lang NP *Clin Oral Implants Res.* 2012 Oct;23 Suppl 6:147-56.
- 10 Cardaropoli D, et al. *Int J Periodontics Restorative Dent.* 2012 Aug;32(4):421-30.
- 11 Cardaropoli D, et al. *Int J Periodontics Restorative Dent.* 2014 Mar-Apr;34(2):211-7.
- 12 Morjaria KR, et al. *Clin Implant Dent Relat Res.* 2014 Feb;16(1):1-20.
- 13 Horváth A, et al. *Clin Oral Investig.* 2013 Mar;17(2):341-63.
- 14 Vittorini Orgeas G, et al. *Int J Oral Maxillofac Implants.* 2013 Jul-Aug;28(4):1049-61.
- 15 Vignoletti F, et al. *Clin Oral Implants Res.* 2012 Feb;23 Suppl 5:22-38.
- 16 Weng D, et al. *Eur J Oral Implantol.* 2011;4 Suppl:59-66.
- 17 Ackermann KL. *Int J Periodontics Restorative Dent.* 2009 Oct;29(5):489-97.
- 18 Schlee M & Esposito M. *Eur J Oral Implantol.* 2009 Autumn;2(3):209-17.
- 19 Pictures by courtesy of Dr. Fernán López
- 20 Picture by courtesy of Prof. Anton Sculean
- 21 Pictures by courtesy of Dr. Juanjo Iturralde Jr.

# RIDGE PRESERVATION WITH GEISTLICH BIOMATERIALS

The use of a biofunctional material such as Geistlich Bio-Oss® is crucial to the long-term successful outcome of extraction socket treatment. After tooth extraction, the slowly resorbing bone matrix Geistlich Bio-Oss® / Geistlich Bio-Oss® Collagen

preserves the ridge volume over time and thus makes a major contribution towards the success of Ridge Preservation<sup>1-3</sup> or ridge contouring at a later time point (e.g. for early implant placement after spontaneous healing)<sup>4,5</sup>

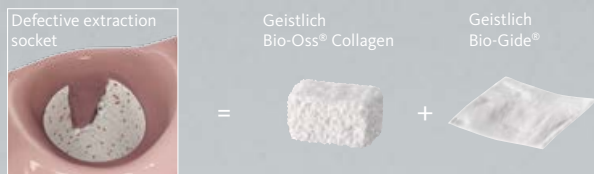

## Clinical benefits of Ridge Preservation with Geistlich Bio-Oss®

Clinical studies indicate that Ridge Preservation using Geistlich Bio-Oss® allows for:

- > stable crest heights in sites with thin buccal bone walls<sup>6</sup>
- > reduced horizontal bone loss in immediate implantation<sup>7</sup>
- > increased mineralized tissue portion in the socket<sup>8</sup>
- > preserved ridge volume under pontics<sup>9</sup>

## Not all Bone Substitutes are the same – Take a closer look!

In recent controlled clinical trials, Geistlich Bio-Oss® showed:

— better ridge preservation than fast resorbing  $\beta$ -TCP<sup>1</sup>

— better ridge preservation than synthetic hydroxyapatite or gelatine sponge<sup>10</sup>

— more mineralized tissue in sockets than allografts<sup>11</sup>

## Geistlich Bio-Gide® – more new bone<sup>12</sup>

Due to its bilayer structure, the Geistlich Bio-Gide® membrane not only prevents ingrowth of soft tissue, but also acts as a guide for the appropriate early blood vessel<sup>13</sup> development and new bone formation<sup>12</sup>.

### Geistlich Bio-Gide®:

- › allows for uneventful wound healing in an open healing approach<sup>2,3</sup>
- › provides for more new bone formation when combined with Geistlich Bio-Oss® vs. Geistlich Bio-Oss® without membrane<sup>12</sup>

### References

- 1 Jung RE, et al. J Clin Periodontol. 2013 Jan;40(1):90-8.
  - 2 Cardaropoli D, et al. Int J Periodontics Restorative Dent. 2012 Aug;32(4):421-30.
  - 3 Cardaropoli D, et al. Int J Periodontics Restorative Dent. 2014 Mar-Apr;34(2):211-7.
  - 4 Buser D, et al. J Dent Res. 2013 Dec;92(12 Suppl):1765-82S.
  - 5 Jensen SS, et al. J Periodontol. 2014 Nov;85(11):1549-56.
  - 6 Nevins M, et al. Int J Periodontics Restorative Dent. 2006 Feb;26(1):19-29.
  - 7 Chen ST, et al. Clin Oral Implants Res. 2007 Oct;18(5):552-62.
  - 8 Lindhe J, et al. Clin Oral Implants Res. 2014 Jul;25(7):786-90.
  - 9 Schlee M & Esposito M. Clin Oral Implants Res. 2014 Jul;25(7):786-90.
  - 10 Shakibaie-M B. Int J Periodontics Restorative Dent. 2013 Mar-Apr;33(2):223-8.
  - 11 Lee DW, et al. Int J Oral Maxillofac Implants. 2009 Jul-Aug;24(4):609-15.
  - 12 Perelman-Karmon M, et al. Int J Periodontics Restorative Dent. 2012 Aug;32(4):459-65.
  - 13 Rothamel D, et al. Clin. Oral Implants Res. 2005;16:369-378.
  - 14 Geistlich Mucograft® Seal Advisory Board Meeting Report 2013. Data on file, Geistlich Pharma AG, Wolhusen, Switzerland.
  - 15 Thoma DS, et al. J Clin Periodontol. 2012 Feb;39(2):157-65.
- \* The definition of an intact extraction socket varies among experts and includes buccal bone defects of 0 to 50%.

## Seal the socket

The collagen matrix of Geistlich Mucograft® Seal specially designed for soft-tissue regeneration is recommended to be used in combination with Geistlich Bio-Oss® Collagen after tooth extraction, when the alveolar buccal walls are preserved.<sup>14</sup>



### Clinical data demonstrates that Geistlich Mucograft® Seal:

- › may enhance early wound healing<sup>15</sup>
- › in combination with Bio-Oss® Collagen significantly reduces the bone loss when compared to spontaneous healing<sup>1</sup>
- › offers flexibility in the therapy concepts: from early implantation 8–10 weeks after extraction through to late implantation or bridge restoration<sup>14</sup>

## Conclusion

- › + 93% ridge width maintained with Geistlich Bio-Oss® Collagen and Geistlich Bio-Gide®<sup>2,3</sup>
- › + 83% ridge width maintained with Geistlich Bio-Oss® Collagen and Geistlich Mucograft® Seal<sup>1</sup>

In the following pages you will find a collection of documented clinical cases showing a great variety of treatment concepts with different Biomaterials.

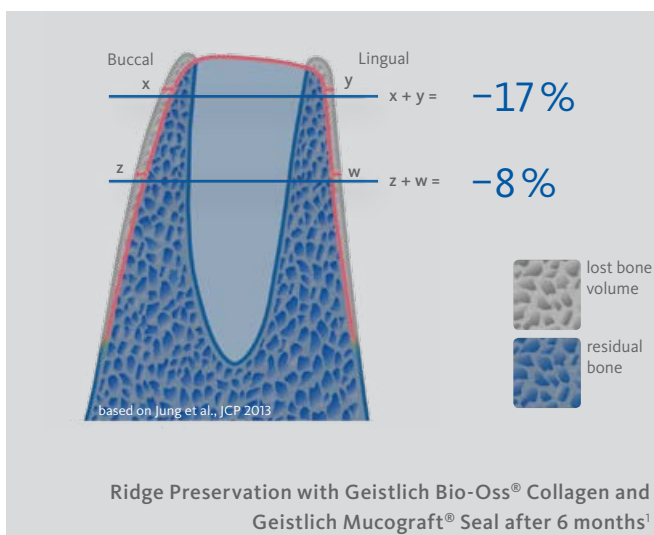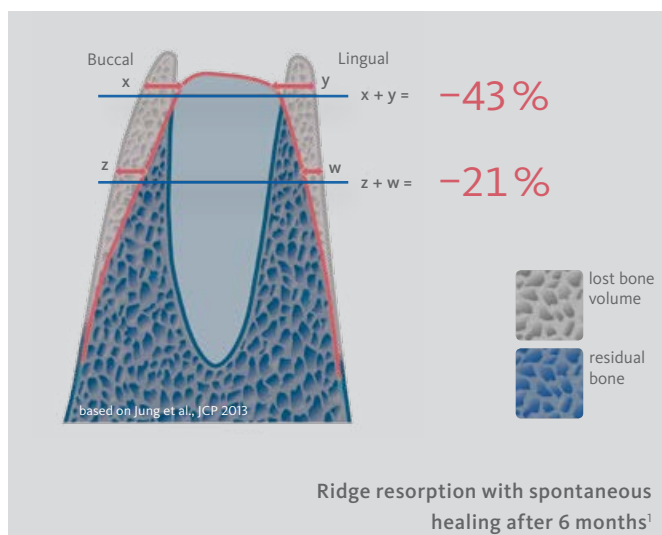

# IMMEDIATE IMPLANT PLACEMENT WITH MINOR BONY DEFECT

## The patient's risk profile

Aesthetic risk factors	Low risk	Medium risk	High risk
Patient's health	Intact immune system (non-smoker)	Light smoker	Impaired immune system (heavy smoker)
Patient's aesthetic requirements	Low	Medium	High
Height of the smile line	Low	Medium	High
Gingival biotype	Thick "low scalloped"	Medium "medium scalloped"	Thin "high scalloped"
Shape of dental crowns	Rectangular		Triangular
Infections at implantation site	None	Chronic	Acute
Bone height at adjacent tooth	≤ 5 mm from contact point	5.5–6.5 mm from contact point	≥ 7 mm from contact point
Restorative status of adjacent tooth	Intact		Restored
Width of tooth gap	1 tooth (≥ 7 mm)	1 tooth (<7mm)	2 teeth or more
Soft-tissue anatomy	Intact		Defective
Bone anatomy of the alveolar ridge	No defect	Horizontal defect	Vertical defect

## Quintessence

### Objectives

- › Immediate implant placement maintaining a good aesthetic contour
- › Minimally invasive procedure
- › The materials are unlimited and easy to use

### Conclusions

- › One-stage treatment of hard and soft tissues
- › The three combined Geistlich Biomaterials heal uneventfully and preserve the alveolar bone volume
- › Long-term good aesthetic outcome in front teeth with short treatment time.








“Geistlich Biomaterials are innovative materials and products and have many years of experience.”

Case documentation

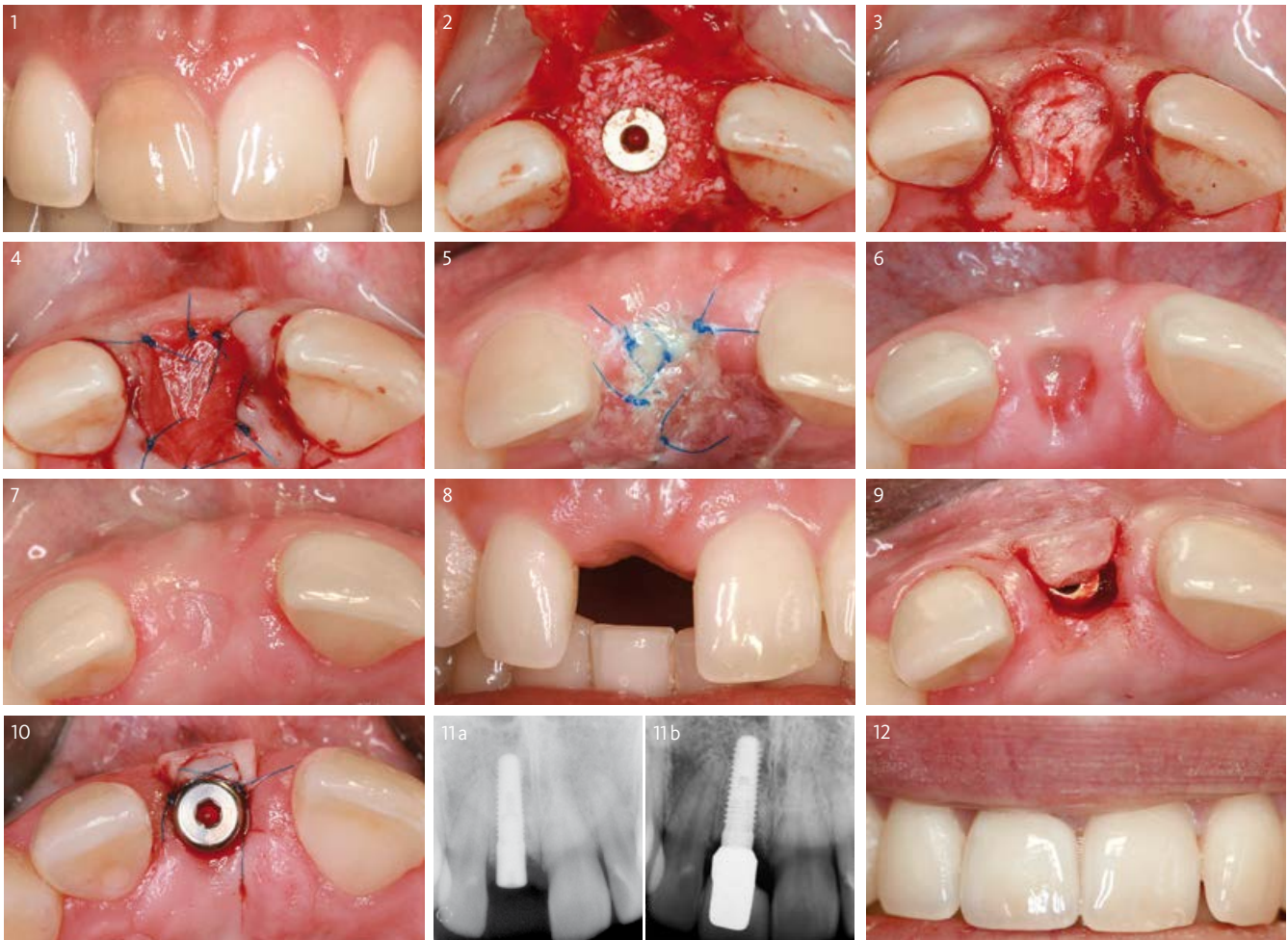

- 1 Clinical situation before extraction of tooth 11.
- 2 Immediate implant placement after flap elevation.
- 3 Implant and Geistlich Bio-Oss® covered with Geistlich Bio-Gide®.
- 4 Geistlich Mucograft® sutured with single sutures on top of augmented area.
- 5 Clinical situation 2 days post-op.
- 6 Clinical situation 3 weeks after surgery.
- 7 Clinical situation 6 weeks after surgery (occlusal).
- 8 Clinical situation 6 weeks after surgery (buccal).
- 9 Minimally invasive re-entry: roll flap preparation after 12 weeks.
- 10 The flap is rolled buccally to increase the thickness of the soft tissues in the buccal area and the abutment is connected.
- 11 a) X-ray after roll flap, 3 months after extraction. b) X-ray at 1-year follow-up.
- 12 Clinical situation 1 year after extraction.

Material selection



Geistlich Bio-Oss® small granules (0.25–1 mm)  
 Geistlich Bio-Gide® (25 × 25 mm)  
 Geistlich Mucograft® (15 × 20 mm; punch 8 mm diameter)

# IMMEDIATE IMPLANT PLACEMENT WITH FILL THE GAP

## The patient's risk profile

Aesthetic risk factors	Low risk	Medium risk	High risk
Patient's health	Intact immune system (non-smoker)	Light smoker	Impaired immune system (heavy smoker)
Patient's aesthetic requirements	Low	Medium	High
Height of the smile line	Low	Medium	High
Gingival biotype	Thick "low scalloped"	Medium "medium scalloped"	Thin "high scalloped"
Shape of dental crowns	Rectangular		Triangular
Infections at implantation site	None	Chronic	Acute
Bone height at adjacent tooth	≤ 5 mm from contact point	5.5–6.5 mm from contact point	≥ 7 mm from contact point
Restorative status of adjacent tooth	Intact		Restored
Width of tooth gap	1 tooth (≥ 7 mm)	1 tooth (<7mm)	2 teeth or more
Soft-tissue anatomy	Intact		Defective
Bone anatomy of the alveolar ridge	No defect	Horizontal defect	Vertical defect

## Quintessence

### Objectives

- › Immediate implant placement in order to reduce the treatment period for the patient
- › Preservation of the vestibular bone volume
- › Preservation of the gingival architecture

### Conclusions

- › The technique minimises the treatment time
- › The treatment maintains the archetype of the soft and hard tissues





Case documentation



- 1 The patient presents with a fractured central incisor. The biotype is rather thin with scalloped marginal gingiva.
- 2 a) X-ray of the fractured tooth. b) Analysis of the bony situation through CBCT allows planning of Type 1 implant placement.
- 3 The gap from implant to the buccal bone is filled with Geistlich Bio-Oss®. A connective tissue graft is placed between the mucosa and the buccal bone.
- 4 The implant (NobelActive™) is positioned optimally, with a more palatal vestibular orientation. The provisional abutment is placed.
- 5 An ideal emergence profile is effected. The provisional crown allows maintenance of the papillae.
- 6 The provisional prosthesis is placed and left out of occlusion.
- 7 Clinical situation 8 days post-operative. The healing occurs uneventfully.
- 8 Situation 4 months after extraction, prior to finalising the prosthetic restoration.
- 9 The natural profile of the soft tissues has been preserved.
- 10 An individual impression post is created for precise transfer of the emergence profile to the lab.
- 11 The final crown is made directly over a zirconia abutment (Procera®).
- 12 Vestibular view of the final restoration 12 months after tooth extraction. Note the perfect alignment of the neck of the teeth and ideal position of the papillae in relation to the contact points.

Material selection



Geistlich Bio-Oss® small granules (0.25–1 mm)

# EARLY IMPLANT PLACEMENT WITH GBR AFTER 8 WEEKS OF SPONTANEOUS HEALING

## The patient's risk profile

Aesthetic risk factors	Low risk	Medium risk	High risk
Patient's health	Intact immune system (non-smoker)	Light smoker	Impaired immune system (heavy smoker)
Patient's aesthetic requirements	Low	Medium	High
Height of the smile line	Low	Medium	High
Gingival biotype	Thick "low scalloped"	Medium "medium scalloped"	Thin "high scalloped"
Shape of dental crowns	Rectangular		Triangular
Infections at implantation site	None	Chronic	Acute
Bone height at adjacent tooth	≤ 5 mm from contact point	5.5–6.5 mm from contact point	≥ 7 mm from contact point
Restorative status of adjacent tooth	Intact		Restored
Width of tooth gap	1 tooth (≥ 7 mm)	1 tooth (< 7mm)	2 teeth or more
Soft-tissue anatomy	Intact		Defective
Bone anatomy of the alveolar ridge	No defect	Horizontal defect	Vertical defect

## Quintessence

### Objectives

- › Pleasing aesthetic outcome
- › Long-term stable bone and soft-tissue situation in the aesthetic region


### Conclusions

- › The low substitution rate of Geistlich Bio-Oss® helps to maintain the volume of the alveolar ridge over time, which is crucial for the long-term aesthetic outcome.
- › Minimal marginal bone loss and low risk of mucosal recession.





Case documentation

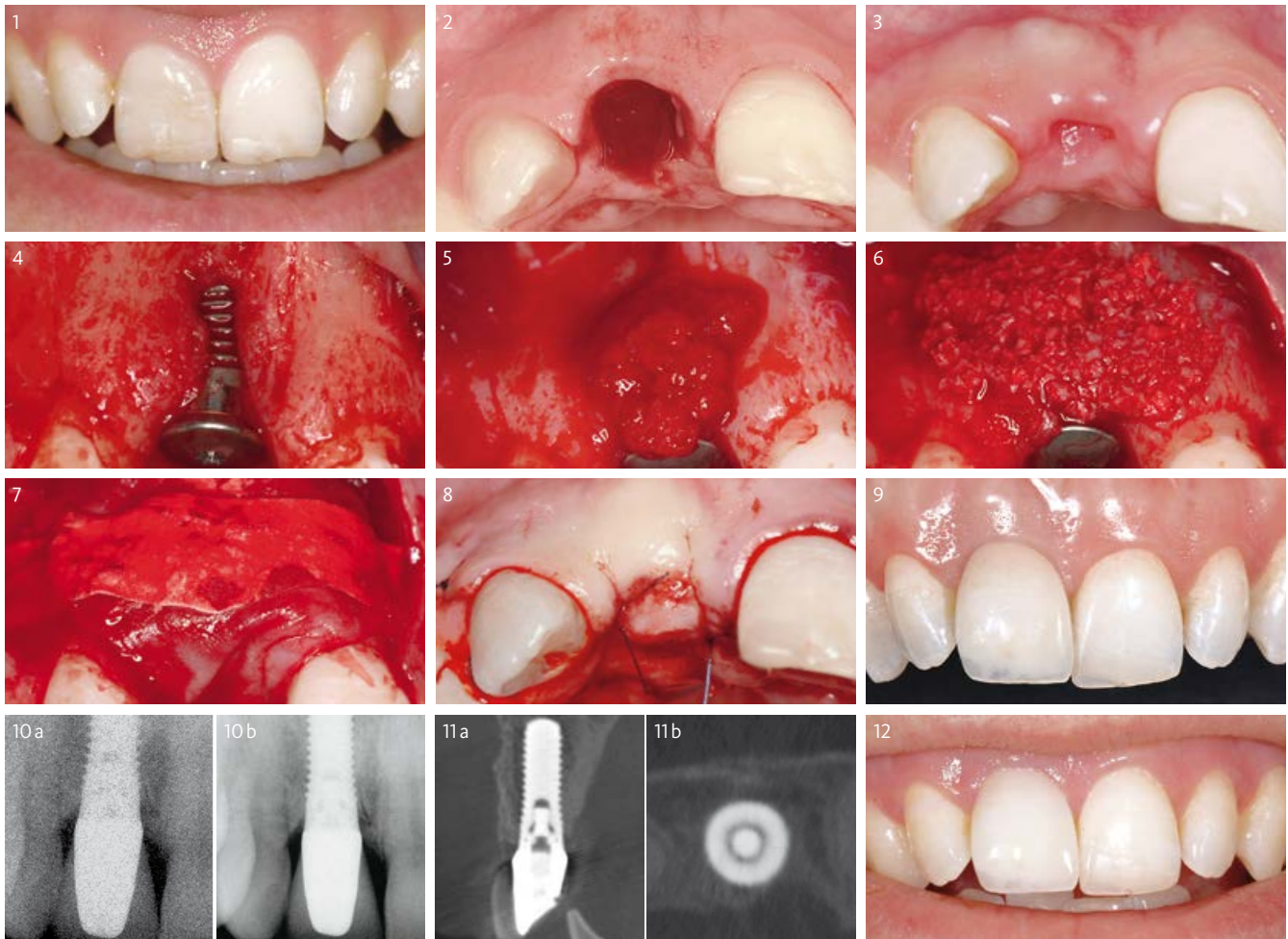

- 1 Clinical findings in the initial examination. The patient exhibits a high smile line and reports an accident several years ago, which affected tooth 11.
- 2 The extraction socket and the soft tissue are allowed to heal for 4–8 weeks after debridement of the inflammatory tissue.
- 3 Within 4–8 weeks of soft tissue healing, no reduction is visible in the crest width in the approximal region of the socket.
- 4 Special attention is paid to correct prosthetic positioning of the implant in all three dimensions with good primary stability.
- 5 The defect is covered with locally harvested autogenous bone chips to promote new bone formation as quickly as possible.
- 6 The bone volume is further optimised by local augmentation using Geistlich Bio-Oss® granules.
- 7 Geistlich Bio-Gide® is applied in two layers to act as a temporary barrier and as a stabiliser for the graft.
- 8 Following the release of the flap by means of mucoperiosteal incisions, a tension-free primary wound closure is achieved. Provisional implant prosthesis starts after 8 weeks.
- 9 The 7.5-year follow-up shows a stable aesthetic outcome.
- 10 X-rays a) at 1 year: implant optimally integrated in the bone; b) at 4 years: absolutely stable peri-implant bony conditions.
- 11 CBCT findings at 7.5 years a) section showing a completely intact facial wall; b) 3-dimensionally correctly placed implant.
- 12 The long-term aesthetic result is excellent.

Material selection



Geistlich Bio-Oss® small granules (0.25–1 mm)  
 Geistlich Bio-Gide® (25 × 25 mm)

# SPONTANEOUS HEALING FOR CANTILEVER IMPLANT BRIDGE

## The patient's risk profile

Aesthetic risk factors	Low risk	Medium risk	High risk
Patient's health	Intact immune system (non-smoker)	Light smoker	Impaired immune system (heavy smoker)
Patient's aesthetic requirements	Low	Medium	High
Height of the smile line	Low	Medium	High
Gingival biotype	Thick "low scalloped"	Medium "medium scalloped"	Thin "high scalloped"
Shape of dental crowns	Rectangular		Triangular
Infections at implantation site	None	Chronic	Acute
Bone height at adjacent tooth	≤ 5 mm from contact point	5.5–6.5 mm from contact point	≥ 7 mm from contact point
Restorative status of adjacent tooth	Intact		Restored
Width of tooth gap	1 tooth (≥ 7 mm)	1 tooth (< 7mm)	2 teeth or more
Soft-tissue anatomy	Intact		Defective
Bone anatomy of the alveolar ridge	No defect	Horizontal defect	Vertical defect

## Quintessence

### Objectives

- › Prosthetic restoration of 2 side-by-side sockets in the anterior area
- › Ridge Preservation for cantilever implant bridge

### Conclusions


- › Early implant placement is suitable for 2 side-by-side sockets
- › The collapse of the tissues during the 6-week healing period can be compensated with a GBR contouring with Geistlich Bio-Oss® and Geistlich Bio-Gide®.





“Early implantation with simultaneous contour augmentation is predictable in challenging cases in the aesthetic zone.”

Case documentation

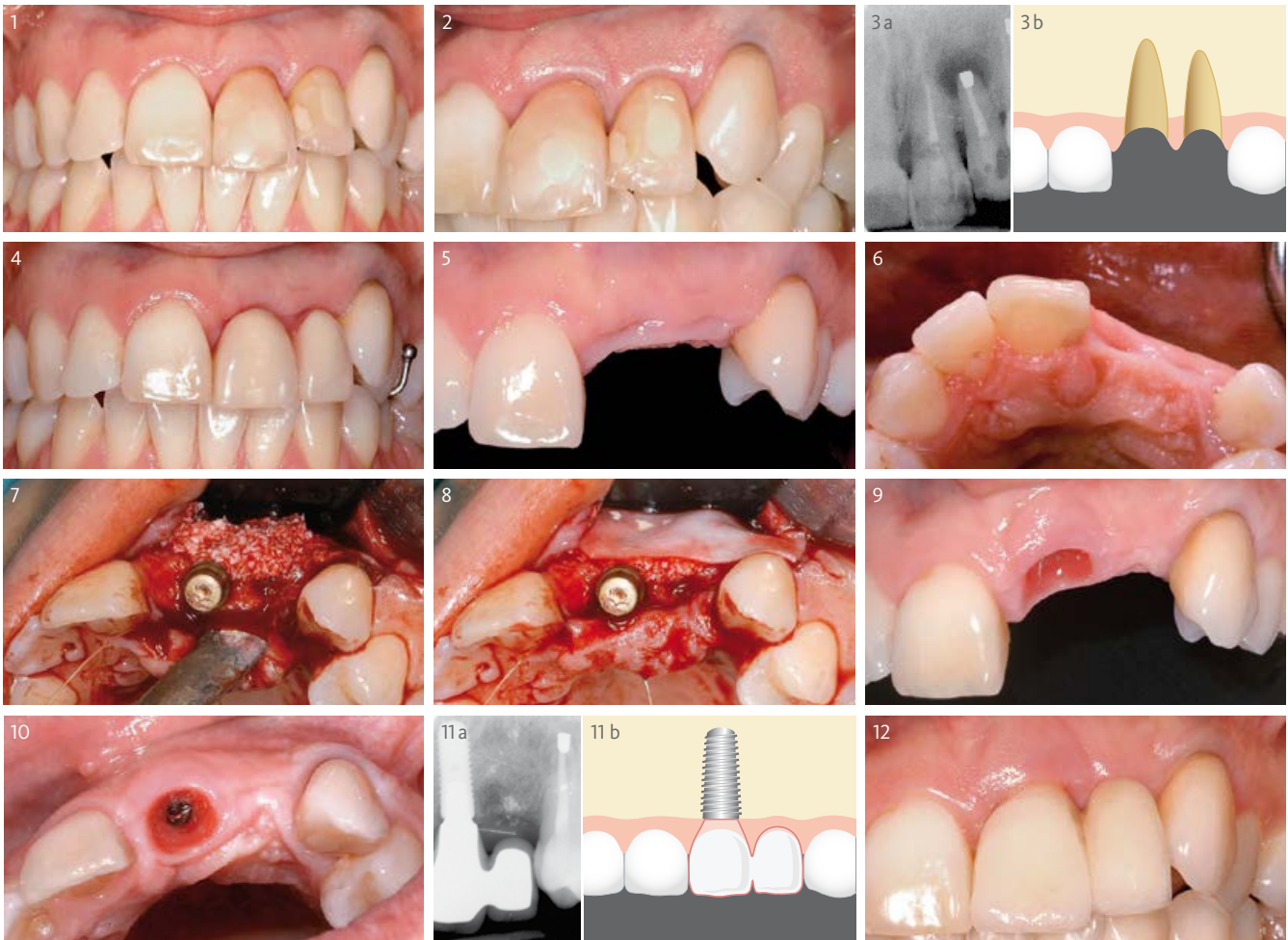

- 1 Initial situation before extraction of 21 and 22.
- 2 Clinical close-up of the pre-operative site prior to extraction of the teeth.
- 3 a) Radiographic findings of the pre-operative site. Note the apical bone resorption at 22 and internal root resorption of tooth 21. b) Scheme of the 2 side-by-side sockets.
- 4 Teeth 21 and 22 are extracted and left heal spontaneously under a provisional restoration.
- 5 Buccal view after 6 weeks of spontaneous healing. Immediately before re-entry. Note the flattening of the ridge anticipating a horizontal defect.
- 6 Occlusal view 6 weeks post-extraction. The soft tissues are healed.
- 7 After flap elevation and implant placement, the resorption of the alveolar bone is compensated with Geistlich Bio-Oss®.
- 8 Geistlich Bio-Gide® is placed over the treated area to stabilise the graft and to obtain the desired contour augmentation.
- 9 Healing of the treated site 18 weeks post-extraction.
- 10 Occlusal view after 18 weeks. Transmucosal healing took place with conditioning of the soft tissues with the provisional crown. The recession on tooth 23 has been covered with a coronally advanced flap and a connective tissue graft.
- 11 a) X-ray of the final prosthetic restoration. b) Schematic representation of the cantilever implant bridge.
- 12 Final situation with the cantilever implant bridge in place 5.5 months after tooth extraction.

Material selection



Geistlich Bio-Oss® small granules (0.25–1 mm)  
 Geistlich Bio-Gide® (25 × 25 mm)

# EARLY IMPLANT PLACEMENT IN EXTRACTION SOCKET WITH PRESERVED BONE WALLS

## The patient's risk profile

Aesthetic risk factors	Low risk	Medium risk	High risk
Patient's health	Intact immune system (non-smoker)	Light smoker	Impaired immune system (heavy smoker)
Patient's aesthetic requirements	Low	Medium	High
Height of the smile line	Low	Medium	High
Gingival biotype	Thick "low scalloped"	Medium "medium scalloped"	Thin "high scalloped"
Shape of dental crowns	Rectangular		Triangular
Infections at implantation site	None	Chronic	Acute
Bone height at adjacent tooth	≤ 5 mm from contact point	5.5–6.5 mm from contact point	≥ 7 mm from contact point
Restorative status of adjacent tooth	Intact		Restored
Width of tooth gap	1 tooth (≥ 7 mm)	1 tooth (<7mm)	2 teeth or more
Soft-tissue anatomy	Intact		Defective
Bone anatomy of the alveolar ridge	No defect	Horizontal defect	Vertical defect

## Quintessence

### Objectives

- › Compensation of the bone resorption through Ridge Preservation
- › Provide the patient with a final restoration in a relatively short time period of time

### Conclusions


- › Almost complete maintenance of the ridge volume is achieved
- › After 8–10 weeks, the soft tissue has a quality and maturity that is adequate for early implant placement.







Case documentation

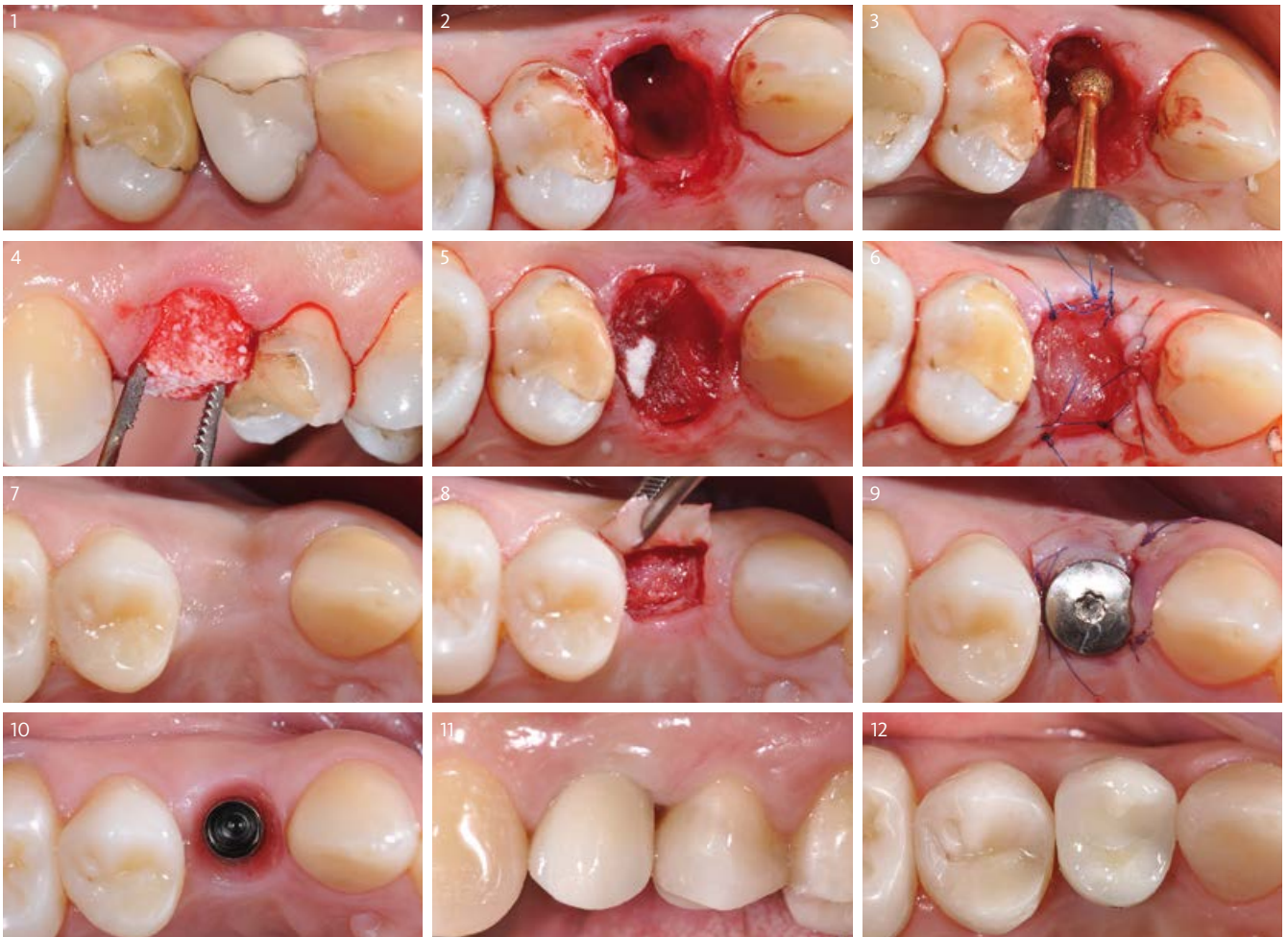

- 1 Initial situation before extraction of tooth 14.
- 2 No buccal bone defect is detected after tooth extraction.
- 3 Extraction socket with de-epithelialised wound margins.
- 4 Extraction socket filled with Geistlich Bio-Oss® Collagen.
- 5 The extraction socket is sealed with Geistlich Mucograft® Seal.
- 6 Geistlich Mucograft® Seal sutured with single interrupted sutures.
- 7 Pre-op clinical situation 10 weeks after extraction (prior to implant placement).
- 8 Preparation of a minimally invasive flap.
- 9 Implant placement with a minimally invasive roll flap technique to improve soft-tissue thickness at the buccal aspect.
- 10 Clinical situation of the soft tissues 4 months after implant placement.
- 11 Final restoration 7 months after tooth extraction (buccal).
- 12 Final restoration 7 months after tooth extraction (occlusal).

Material selection



Geistlich Bio-Oss® Collagen (100 mg)  
 Geistlich Mucograft® Seal (8 mm diameter)

# RIDGE PRESERVATION IN SOCKET WITH PRESERVED BUCCAL BONE WALL

## The patient's risk profile

Aesthetic risk factors	Low risk	Medium risk	High risk
Patient's health	Intact immune system (non-smoker)	Light smoker	Impaired immune system (heavy smoker)
Patient's aesthetic requirements	Low	Medium	High
Height of the smile line	Low	Medium	High
Gingival biotype	Thick "low scalloped"	Medium "medium scalloped"	Thin "high scalloped"
Shape of dental crowns	Rectangular		Triangular
Infections at implantation site	None	Chronic	Acute
Bone height at adjacent tooth	≤ 5 mm from contact point	5.5–6.5 mm from contact point	≥ 7 mm from contact point
Restorative status of adjacent tooth	Intact		Restored
Width of tooth gap	1 tooth (≥ 7 mm)	1 tooth (<7mm)	2 teeth or more
Soft-tissue anatomy	Intact		Defective
Bone anatomy of the alveolar ridge	No defect*	Horizontal defect	Vertical defect

## Quintessence

### Objectives

- › Preservation of hard and soft-tissue volume after tooth extraction.
- › Late implant placement, as it is an extremely reliable procedure, which has been proven repeatedly in the international literature.

### Conclusions

- › Geistlich Bio-Oss® Collagen and Geistlich Mucograft® Seal preserve the ridge for optimal implant placement 5 months post-op.
- › At the central incisor, the buccal soft-tissue thickness is optimised with a connective tissue graft.




\* Buccal bone wall preserved, but more apically with respect to the neighbouring teeth because of a discrepancy on the marginal gingiva level.



“With the chosen Biomaterials, hard and soft-tissue volume are preserved in the front area for late implantation.”

Case documentation

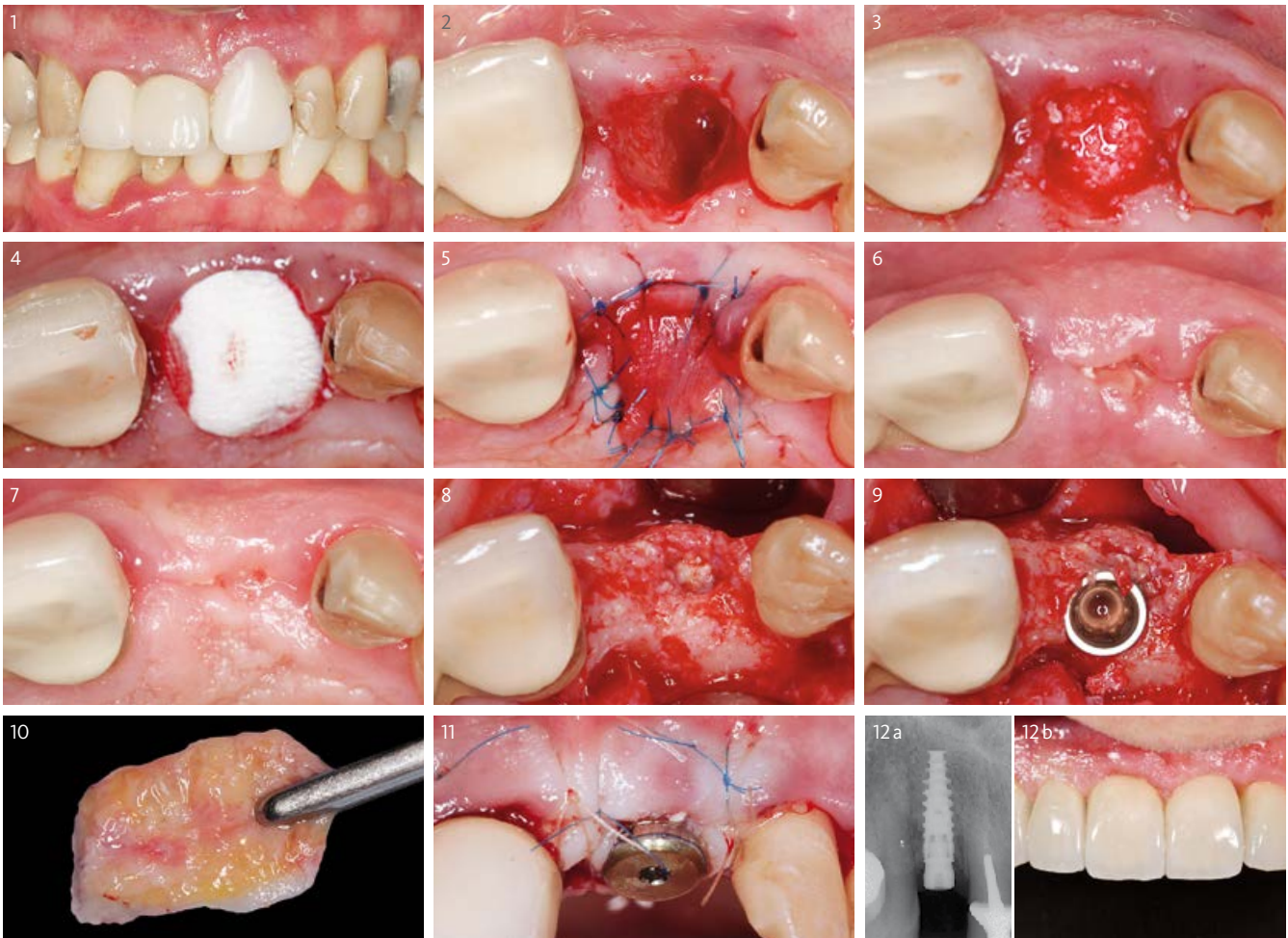

- 1 Tooth 21 is scheduled for extraction due to periodontal problems.
- 2 Meticulous curettage of the socket after atraumatic flap less extraction.
- 3 Filling of the extraction socket with Geistlich Bio-Oss® Collagen up to the palatal bone.
- 4 Geistlich Mucograft® Seal in place: the spongy structure faces towards the bone substitute.
- 5 Geistlich Mucograft® Seal is sutured with single interrupted sutures allowing optimal adaptation between the borders of the soft tissues and the collagen matrix.
- 6 Wound healing at 2 weeks: good healing of the soft tissues with a beautiful pink colour.
- 7 Wound healing at 3 months: complete closure of the socket with mature soft tissues.
- 8 Five months after extraction: good maintenance of the alveolar bone volume.
- 9 Implant placement to replace tooth 21 without additional GBR.
- 10 Connective tissue graft harvested at the left palate.
- 11 The connective tissue graft is placed at the buccal site and the flap is closed with suspension sutures and single interrupted sutures (monofilament 6/0).
- 12 a) X-ray shows the osseointegrated implant 3 months after implant placement. b) Follow-up 28 months after extraction.

Material selection



Geistlich Bio-Oss® Collagen (100 mg)  
 Geistlich Mucograft® Seal (8 mm diameter)

# RIDGE PRESERVATION IN EXTRACTION SOCKET WITH PRESERVED BUCCAL BONE

## The patient's risk profile

Aesthetic risk factors	Low risk	Medium risk	High risk
Patient's health	Intact immune system (non-smoker)	Light smoker	Impaired immune system (heavy smoker)
Patient's aesthetic requirements	Low	Medium	High
Height of the smile line	Low	Medium	High
Gingival biotype	Thick "low scalloped"	Medium "medium scalloped"	Thin "high scalloped"
Shape of dental crowns	Rectangular		Triangular
Infections at implantation site	None	Chronic	Acute
Bone height at adjacent tooth	≤ 5 mm from contact point	5.5–6.5 mm from contact point	≥ 7 mm from contact point
Restorative status of adjacent tooth	Intact		Restored
Width of tooth gap	1 tooth (≥ 7 mm)	1 tooth (< 7mm)	2 teeth or more
Soft-tissue anatomy	Intact		Defective
Bone anatomy of the alveolar ridge	No defect	Horizontal defect*	Vertical defect

## Quintessence

### Objectives

- › Delayed implant placement 4 months after extraction
- › Minimally invasive treatment of the socket

### Conclusions

- › Good/mature/solid bone obtained 4 months after treatment
- › Fast and scar-free soft-tissue regeneration
- › Optimal clinical and aesthetic result for the patient



\* Intact extraction socket, with a minor bony defect up to 50% of the buccal bone wall



“Soft and hard tissues are well preserved without any scarring on the buccal or occlusal aspect.”

### Case documentation



- 1 Situation on the day of tooth extraction.
- 2 Pre-op situation (buccal).
- 3 The sulcus is de-epithelialised using a diamond bur.
- 4 The extraction socket is filled with Geistlich Bio-Oss® Collagen.
- 5 Geistlich Mucograft® Seal in place sutured with single and double interrupted sutures.
- 6 Healing of soft tissues 3 days after tooth extraction.

- 7 Healing of the soft tissues at the time of suture removal 10 days after surgery.
- 8 Tissue healing 9 weeks after tooth extraction.
- 9 Situation after 4 months at the time of implant placement.
- 10 The flap elevation reveals ideal bony situation for implant placement.
- 11 Implant seated.
- 12 Final restoration 11 months after tooth extraction.


### Material selection



Geistlich Bio-Oss® Collagen  
(100 mg)  
Geistlich Mucograft® Seal  
(8 mm diameter)

# EXTRACTION SOCKET TREATMENT OPTIONS

The appropriate type of treatment for the management of extraction sockets is derived from a coherent evaluation of the aesthetic risk factors. In addition to the time of implantation, the attending dentist needs to make a decision regarding regenerative measures directly after tooth extraction. Various procedures are recommended:

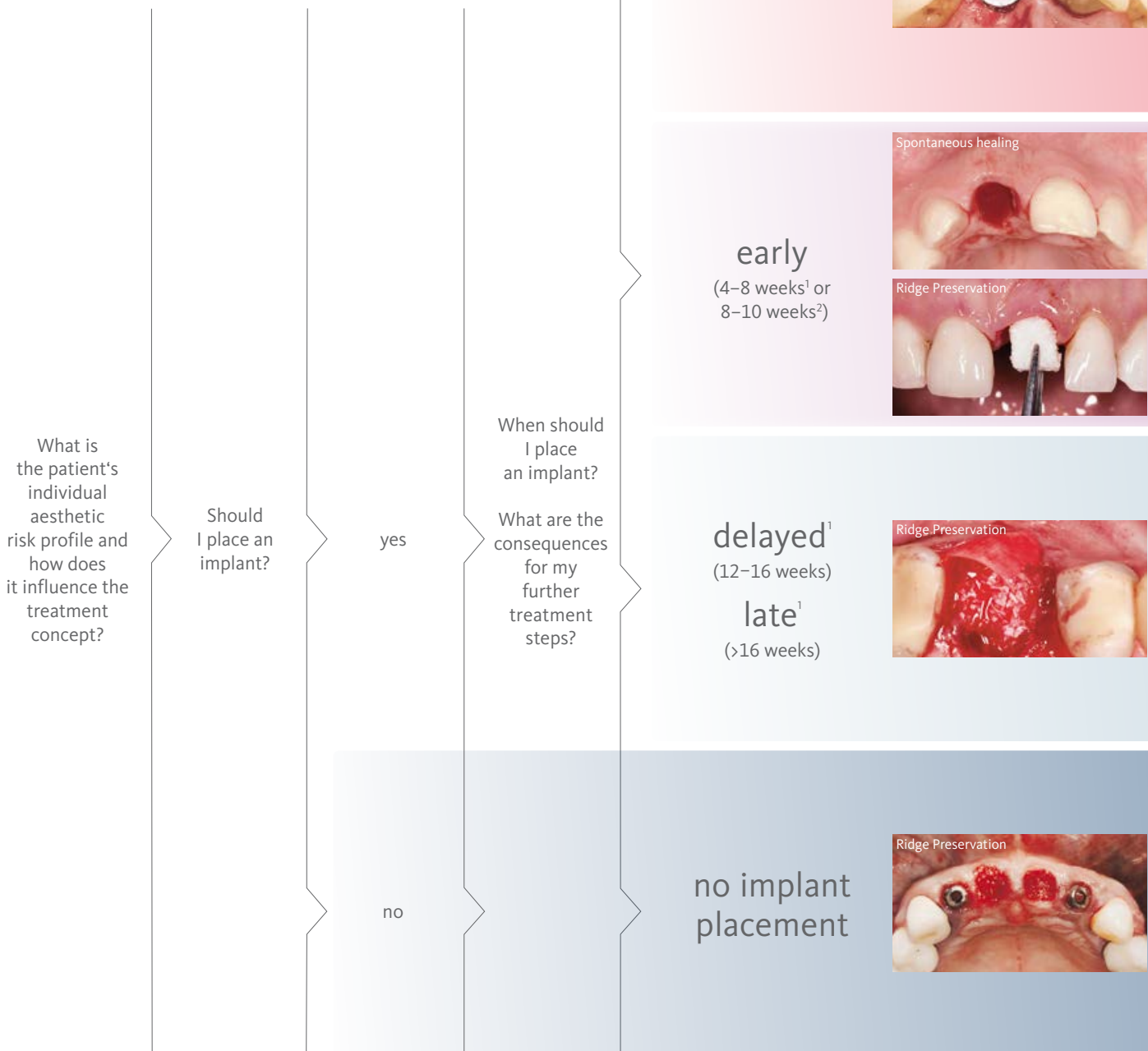

## References

- 1 Hämmerle CH. et al., Int J Oral Maxillofac Implants. 2004;19 Suppl:26-8.
- 2 Geistlich Mucograft® Seal report on the meeting of the Advisory Committee, 2013. Data on file, Geistlich Pharma AG, Wolhusen, Switzerland.

\* The definition of an intact extraction socket varies among experts and includes buccal bone defects of 0 to 50 %.

Timeline

Day 0


8 - 10 weeks

3 months


6 months





Spontaneous healing

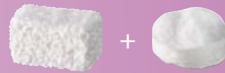

Guided Bone Regeneration at implant placement, if required

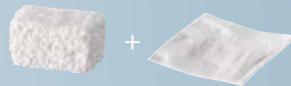

Guided Bone Regeneration at implant placement, if required

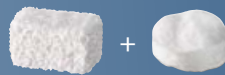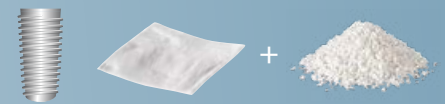

Implant



Geistlich Bio-Oss®



Geistlich Bio-Oss® Collagen



Geistlich Mucograft® Seal



Geistlich Bio-Gide®



Bridge

# RIDGE PRESERVATION IN DEFECT EXTRACTION SOCKETS

## The patient's risk profile

Aesthetic risk factors	Low risk	Medium risk	High risk
Patient's health	Intact immune system (non-smoker)	Light smoker	Impaired immune system (heavy smoker)
Patient's aesthetic requirements	Low	Medium	High
Height of the smile line	Low	Medium	High
Gingival biotype	Thick "low scalloped"	Medium "medium scalloped"	Thin "high scalloped"
Shape of dental crowns	Rectangular		Triangular
Infections at implantation site	None	Chronic	Acute
Bone height at adjacent tooth	≤ 5 mm from contact point	5.5–6.5 mm from contact point	≥ 7 mm from contact point
Restorative status of adjacent tooth	Intact		Restored
Width of tooth gap	1 tooth (≥ 7 mm)	1 tooth (< 7mm)	2 teeth or more
Soft-tissue anatomy	Intact		Defective
Bone anatomy of the alveolar ridge	No defect	Horizontal defect	Vertical defect

## Quintessence

### Objectives

- › Maintain hard and soft-tissue contour in aesthetically demanding region
- › Late implant placement in single tooth gap

### Conclusions

- › Severe ridge resorption was prevented with Geistlich Biomaterials
- › A long-term pleasant outcome was achieved with additional contouring with Geistlich Biomaterials and a connective tissue graft at implant placement







“Whenever possible we prefer to preserve rather than to rebuild the bone later, specially in the front teeth.”

Case documentation



- 1 Initial situation before removal of tooth 21.
- 2 Inspection of the extraction socket with the periodontal probe shows a buccal bony defect.
- 3 Geistlich Bio-Gide® is placed buccally on the inner alveolar wall, slightly protruding the crestal bone. Geistlich Bio-Oss® Collagen fills the socket up to the crestal bone level.
- 4 Geistlich Bio-Oss® (small granules) are packed on top of Geistlich Bio-Oss® Collagen up to soft-tissue level.
- 5 The collagen membrane is folded over the filled socket, adapted under the palatal sulcus, fixed with vertical mattress sutures and heals by secondary intention.
- 6 Uneventful healing situation 3 days post-extraction.
- 7 Clinical situation 1 week after tooth extraction.
- 8 Situation after site-conditioning of the soft tissues 4 months post-extraction.
- 9 Flap elevation and implant placement reveal a fenestration 4 months after tooth extraction.
- 10 The ridge is contoured with a GBR (Geistlich Bio-Oss® and Geistlich Bio-Gide®) and a connective tissue graft on the buccal-crestal area.
- 11 The flap is closed over the graft.
- 12 Loading of the implant with the final restoration 7 months after implant placement (11 months after extraction).

Material selection



Geistlich Bio-Oss® small granules (0.25–1 mm)  
 Geistlich Bio-Oss® Collagen (100 mg)  
 Geistlich Bio-Gide® (25 × 25 mm)

# RIDGE PRESERVATION IN THE ANTERIOR REGION FOR LATE IMPLANTATION

## The patient's risk profile

Aesthetic risk factors	Low risk	Medium risk	High risk
Patient's health	Intact immune system (non-smoker)	Light smoker	Impaired immune system (heavy smoker)
Patient's aesthetic requirements	Low	Medium	High
Height of the smile line	Low	Medium	High
Gingival biotype	Thick "low scalloped"	Medium "medium scalloped"	Thin "high scalloped"
Shape of dental crowns	Rectangular		Triangular
Infections at implantation site	None	Chronic	Acute
Bone height at adjacent tooth	≤ 5 mm from contact point	5.5–6.5 mm from contact point	≥ 7 mm from contact point
Restorative status of adjacent tooth	Intact		Restored
Width of tooth gap	1 tooth (≥ 7 mm)	1 tooth (< 7 mm)	2 teeth or more
Soft-tissue anatomy	Intact		Defective
Bone anatomy of the alveolar ridge	No defect	Horizontal defect*	Vertical defect

## Quintessence

### Objectives

- › Preservation of hard and soft-tissue volume after extraction in the anterior region for late implant placement.
- › Prevention of extensive guided bone regeneration procedures at implant placement.

### Conclusions

- › Volume of hard and soft tissue can be preserved better with Geistlich Bio-Oss® Collagen and Geistlich Mucograft® Seal than with spontaneous healing.<sup>1</sup>
- › A minimally invasive GBR is performed to contour the ridge at implant placement.




<sup>1</sup> Jung RE, et al. J Clin Periodontol. 2013 Jan;40(1):90–8

\* Intact extraction socket, with a minor bony defect up to 50% of the buccal bone wall



Case documentation

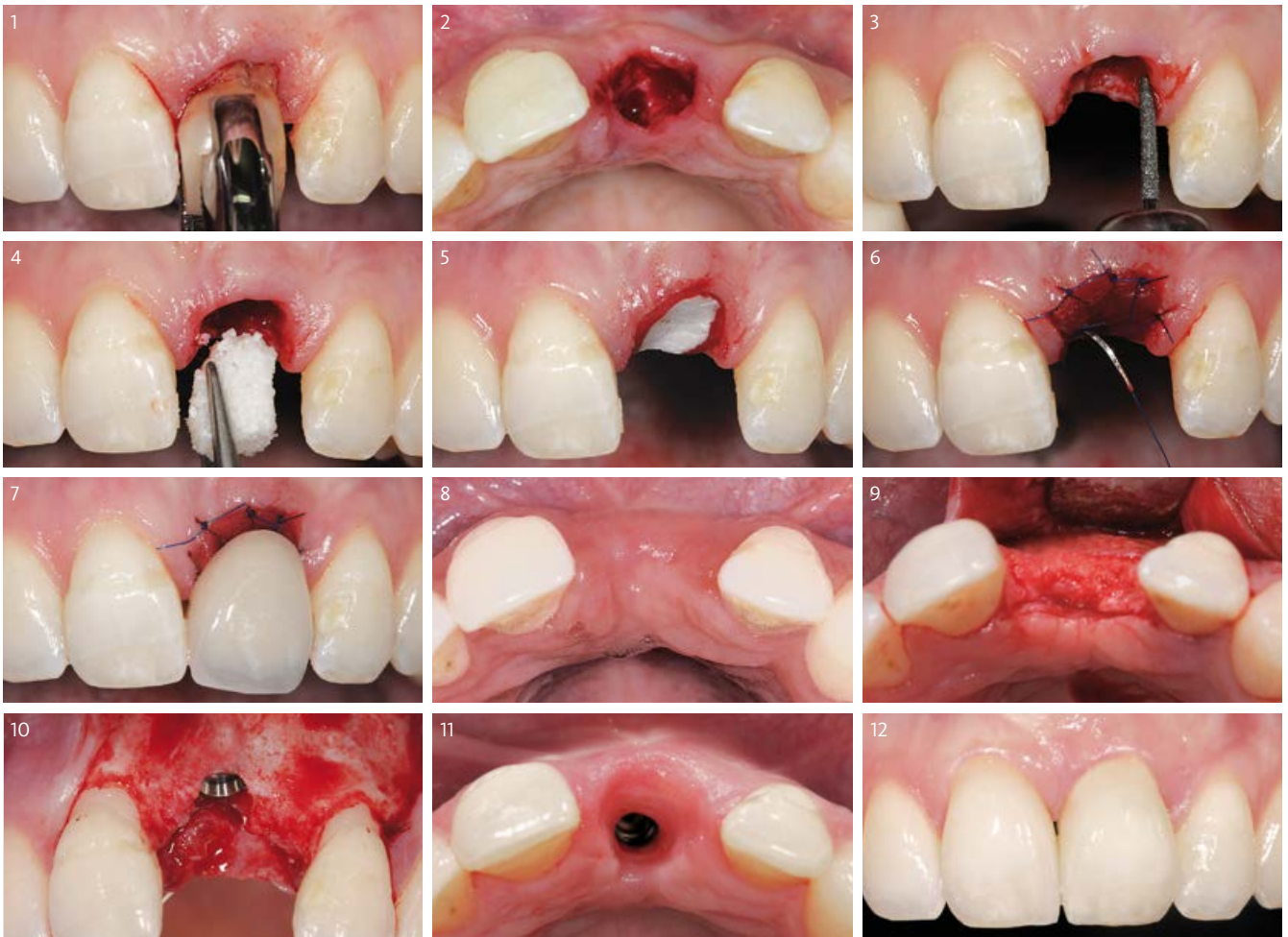

- 1 Extraction of tooth 21 due to a trauma with concomitant external resorptions. Care was taken in preserving the alveolar bone.
- 2 Occlusal view of the socket after tooth extraction. No flaps are raised around the affected area. A slight buccal bone defect was observed.
- 3 The socket is gently curetted for removal of granulation tissue. Subsequently, the wound margins were de-epithelialised with a diamond in a counter-piece with water cooling.
- 4 Filling of the extraction socket with Geistlich Bio-Oss® Collagen to the level of the palatal bone.
- 5 Geistlich Mucograft® is applied dry and adapts perfectly to the wound margins.
- 6 Suturing of the Geistlich Mucograft® with 6-0 single interrupted sutures.
- 7 The tissues are left to heal beneath the provisional, taking care not to apply pressure to the biomaterials.
- 8 Situation 7.5 months after extraction revealing nice soft-tissue situation with a slight dip at the buccal aspect.
- 9 Flap elevation shows the healed bony situation 7.5 months after Ridge Preservation.
- 10 Implant placement in fully mature bone. A small GBR for contouring is performed.
- 11 Excellent emergence profile after 10 months.
- 12 Situation with the final restoration 10 months after tooth extraction.

Material selection



Geistlich Bio-Oss® Collagen (100 mg)  
 Geistlich Mucograft® (15×20 mm punch 8 mm diameter)

# RIDGE PRESERVATION IN THE POSTERIOR REGION FOR LATE IMPLANTATION

## The patient's risk profile

Aesthetic risk factors	Low risk	Medium risk	High risk
Patient's health	Intact immune system (non-smoker)	Light smoker	Impaired immune system (heavy smoker)
Patient's aesthetic requirements	Low	Medium	High
Height of the smile line	Low	Medium	High
Gingival biotype	Thick "low scalloped"	Medium "medium scalloped"	Thin "high scalloped"
Shape of dental crowns	Rectangular		Triangular
Infections at implantation site	None	Chronic	Acute
Bone height at adjacent tooth	≤ 5 mm from contact point	5.5–6.5 mm from contact point	≥ 7 mm from contact point
Restorative status of adjacent tooth	Intact		Restored
Width of tooth gap	1 tooth (≥ 7 mm)	1 tooth (<7 mm)	2 teeth or more
Soft-tissue anatomy	Intact		Defective
Bone anatomy of the alveolar ridge	No defect	Horizontal defect	Vertical defect

## Quintessence

### Objectives

- › Preservation of the ridge contour with minimal invasion
- › Late implant placement

### Conclusions


- › Geistlich Bio-Oss® and Geistlich Mucograft® Seal enable a flapless and effective Ridge Preservation
- › Hard and soft tissues are optimal for implant placement 6 months after Ridge Preservation procedure





## “Geistlich Bio-Oss® and Geistlich Mucograft® Seal enable a flapless and effective Ridge Preservation.”

### Case documentation

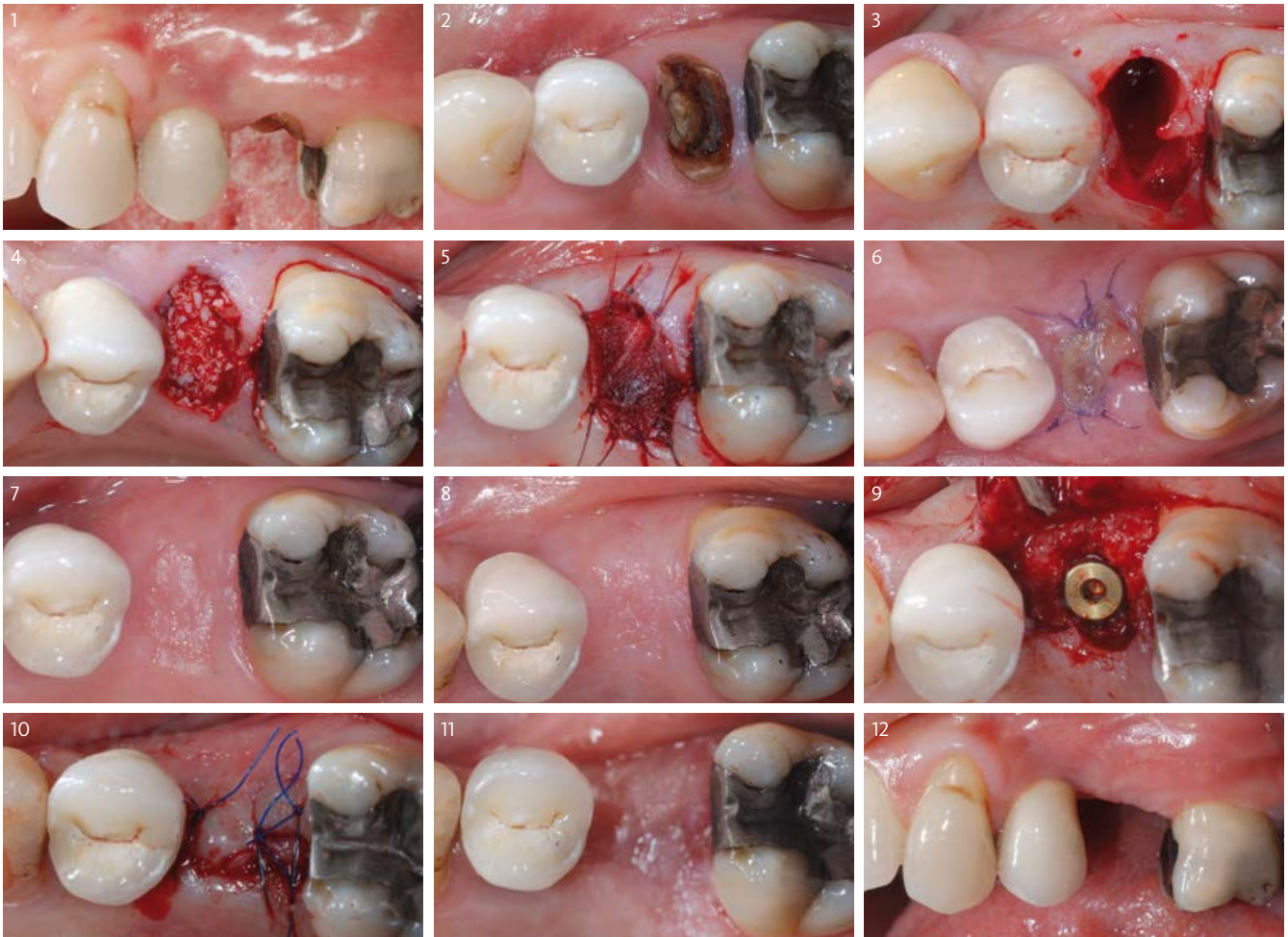

- 1 Clinical appearance before treatment (buccal).
- 2 Clinical appearance before treatment (occlusal).
- 3 Situation after tooth extraction.
- 4 The socket is grafted with Geistlich Bio-Oss® up to the bone level.
- 5 Geistlich Mucograft® Seal is sutured with 8 single interrupted sutures.
- 6 Healing of the soft tissues 1 week after tooth extraction.
- 7 Clinical post-op appearance 8 weeks after extraction.
- 8 Situation 6 months after tooth extraction and before implant placement.
- 9 Minimal flap elevation reveals optimal bony and soft-tissue situation for correct implant placement.
- 10 Closure of the flap for submerged healing.
- 11 Occlusal clinical view 3 weeks after submerged implant placement (6.5 months after extraction).
- 12 Buccal clinical view 6.5 months after extraction.

### Material selection



Geistlich Bio-Oss® small granules (0.25–1 mm)  
 Geistlich Mucograft® Seal (15×20 mm punch 8 mm diameter)

# RIDGE PRESERVATION IN DEFECT EXTRACTION SOCKET

## The patient's risk profile

Aesthetic risk factors	Low risk	Medium risk	High risk
Patient's health	Intact immune system (non-smoker)	Light smoker	Impaired immune system (heavy smoker)
Patient's aesthetic requirements	Low	Medium	High
Height of the smile line	Low	Medium	High
Gingival biotype	Thick "low scalloped"	Medium "medium scalloped"	Thin "high scalloped"
Shape of dental crowns	Rectangular		Triangular
Infections at implantation site	None	Chronic	Acute
Bone height at adjacent tooth	≤ 5 mm from contact point	5.5–6.5 mm from contact point	≥ 7 mm from contact point
Restorative status of adjacent tooth	Intact		Restored
Width of tooth gap	1 tooth (≥ 7 mm)	1 tooth (<7mm)	2 teeth or more
Soft-tissue anatomy	Intact		Defective
Bone anatomy of the alveolar ridge	No defect	Horizontal defect	Vertical defect

## Quintessence


### Objectives

- › Replace a falling hopeless central incisor with high aesthetic demands.

### Conclusions


- › Ridge Preservation techniques are effective in minimising volume loss.



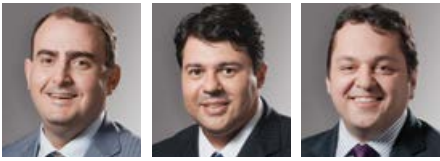

“Geistlich Biomaterials grant us safe choices, effective treatments and predictable outcomes.”

Case documentation

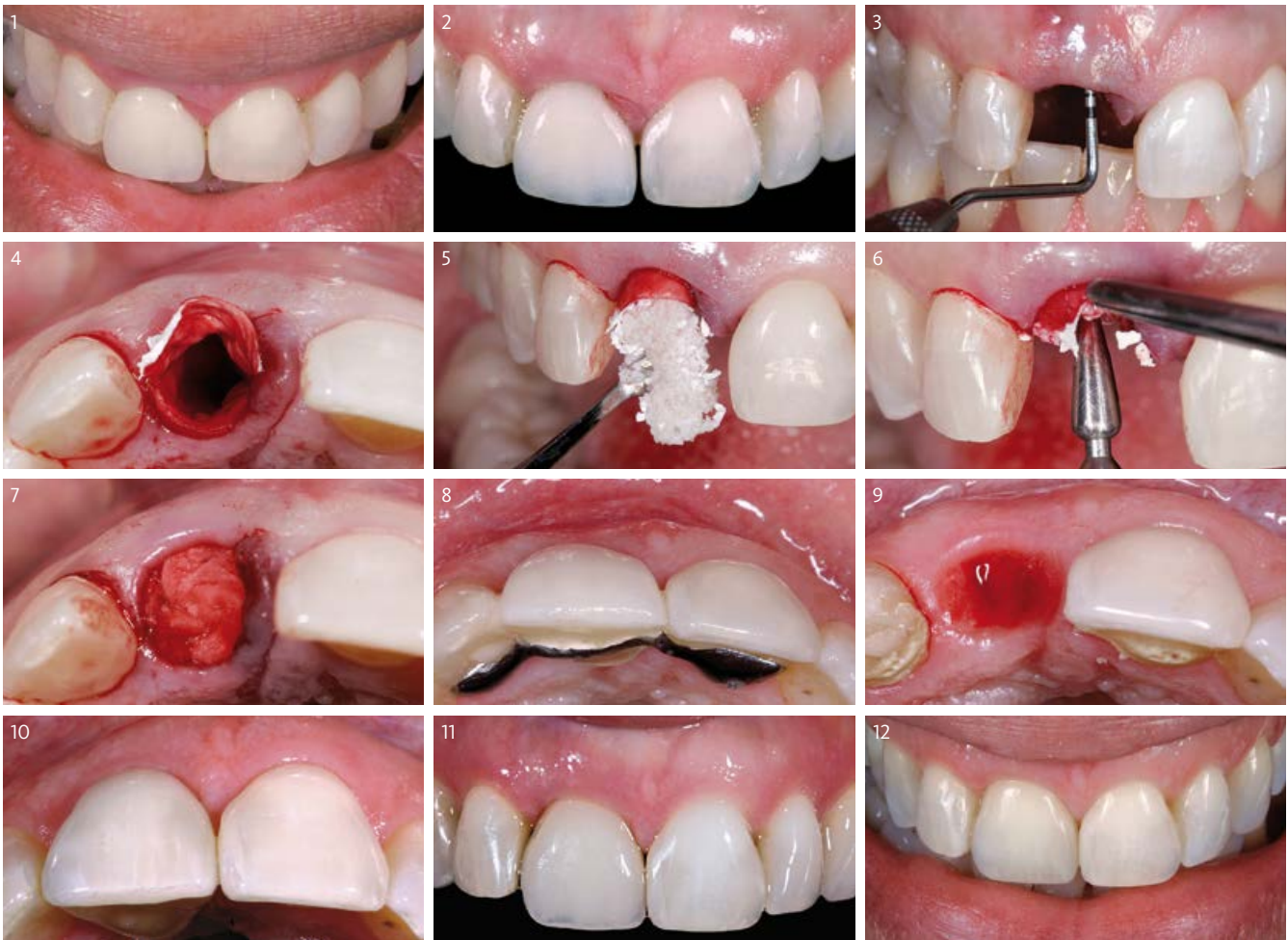

- 1 Initial Smile.
- 2 Initial front view.
- 3 Alveolar inspection after the extraction. Note the presence of deep buccal defect.
- 4 Positioning of Geistlich Bio-Gide® by buccal and palatal.
- 5 Socket filling with Geistlich Bio-Oss®.
- 6 Geistlich Bio-Oss® accommodation .
- 7 Repositioning of Geistlich Bio-Gide® for socket sealing.
- 8 Provisional prosthesis installed.
- 9 Healing aspect after 6 months. Note the volume preservation and tissue contour. Connective tissue graft and late implant placement after 6 months.
- 10 Final rehabilitation with ceramic crown\*. Note the natural contour and emergence profile.
- 11 Final clinical appearance showing harmony in gingival margin positioning.
- 12 Harmonic smile after rehabilitation.

Material selection



Geistlich Bio-Oss® small granules (0.25–1 mm)  
 Geistlich Bio-Gide® (13 × 25 mm)

\* Prosthesis rehabilitation by Prof. Dr. Oswaldo Scopin de Andrade and Luis Alves

# RIDGE PRESERVATION IN DEFECT EXTRACTION SOCKET

## The patient's risk profile

Aesthetic risk factors	Low risk	Medium risk	High risk
Patient's health	Intact immune system (non-smoker)	Light smoker	Impaired immune system (heavy smoker)
Patient's aesthetic requirements	Low	Medium	High
Height of the smile line	Low	Medium	High
Gingival biotype	Thick "low scalloped"	Medium "medium scalloped"	Thin "high scalloped"
Shape of dental crowns	Rectangular		Triangular
Infections at implantation site	None	Chronic	Acute
Bone height at adjacent tooth	≤ 5 mm from contact point	5.5–6.5 mm from contact point	≥ 7 mm from contact point
Restorative status of adjacent tooth	Intact		Restored
Width of tooth gap	1 tooth (≥ 7 mm)	1 tooth (<7mm)	2 teeth or more
Soft-tissue anatomy	Intact		Defective
Bone anatomy of the alveolar ridge	No defect	Horizontal defect	Vertical defect

## Quintessence

### Objectives

- › Prevent tissue collapse in the posterior area due to absence of the buccal bone wall.
- › Avoid a possible sinus lift elevation.

### Conclusions

- › Ridge preservation with Geistlich Biomaterials preserved the alveolar ridge contour.
- › A minimally invasive procedure provided enough ridge width for adequate implant placement and aesthetic outcome.








“Ridge Preservation allows correct 3D implant placement reducing additional surgeries (i.e. sinus lift).”

### Case documentation

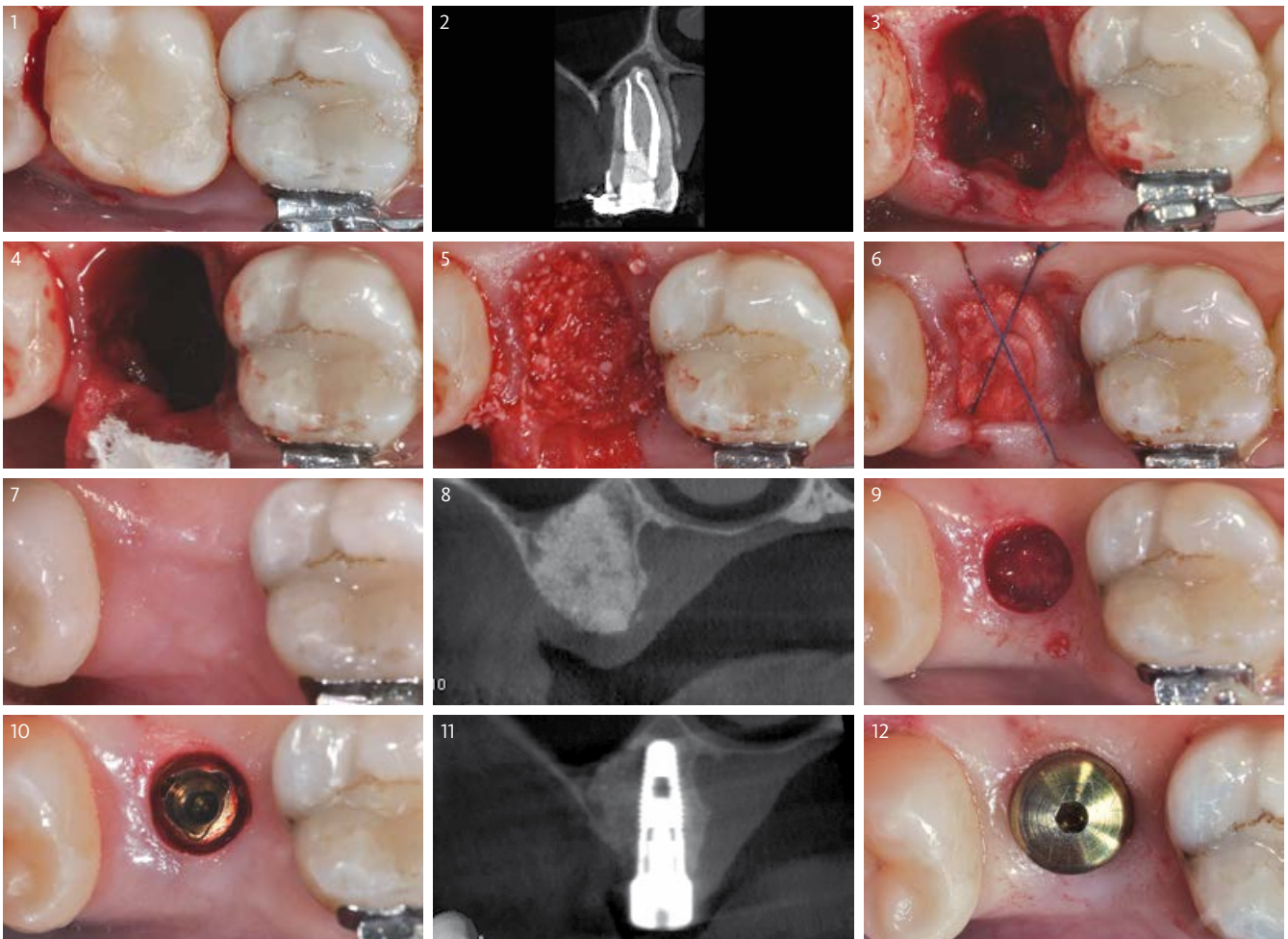

- 1 Compromised upper molar due to longitudinal tooth fracture.
- 2 CBCT upper molar before extraction. Note the absence of the buccal bone wall.
- 3 Socket after tooth extraction.
- 4 Buccal bone wall replaced by Geistlich Bio-Gide®.
- 5 Filling with Geistlich Bio-Oss® (small granules 0.25–1 mm).
- 6 Geistlich Bio-Gide® is sutured with a cross-suture.
- 7 Clinical situation after 6 months of healing.

- 8 CBCT 6 months post-extraction before implant placement.
- 9 Flapless implant installation procedure 6 months after tooth extraction.
- 10 Implant in place 6 months after tooth extraction and Ridge Preservation procedure.
- 11 CBCT immediately after implant placement.
- 12 Abutment connection

### Material selection



Geistlich Bio-Oss® small granules (0.25–1 mm)  
Geistlich Bio-Gide® (25 × 25 mm)

# RIDGE PRESERVATION FOR DELAYED IMPLANT PLACEMENT

## The patient's risk profile

Aesthetic risk factors	Low risk	Medium risk	High risk
Patient's health	Intact immune system (non-smoker)	Light smoker	Impaired immune system (heavy smoker)
Patient's aesthetic requirements	Low	Medium	High
Height of the smile line	Low	Medium	High
Gingival biotype	Thick "low scalloped"	Medium "medium scalloped"	Thin "high scalloped"
Shape of dental crowns	Rectangular		Triangular
Infections at implantation site	None	Chronic	Acute
Bone height at adjacent tooth	≤ 5 mm from contact point	5.5–6.5 mm from contact point	≥ 7 mm from contact point
Restorative status of adjacent tooth	Intact		Restored
Width of tooth gap	1 tooth (≥ 7 mm)	1 tooth (< 7mm)	2 teeth or more
Soft-tissue anatomy	Intact		Defective
Bone anatomy of the alveolar ridge	No defect	Horizontal defect	Vertical defect

## Quintessence

### Objectives

- › Reconstruct alveolar bone with severe vertical loss from chronic periodontitis at the lower left second molar
- › Investigate the clinical and histological result by using Geistlich Combi-Kit Collagen after tooth extraction.

### Conclusions


- › The defect was completely filled with newly-formed hard tissue after 6 months
- › Histomorphometric analysis revealed 45% of the hard tissue area including bone substitute material and 28% of the soft tissue area.





“After 6 months the defect was completely filled with newly-formed hard tissue.”

Case documentation

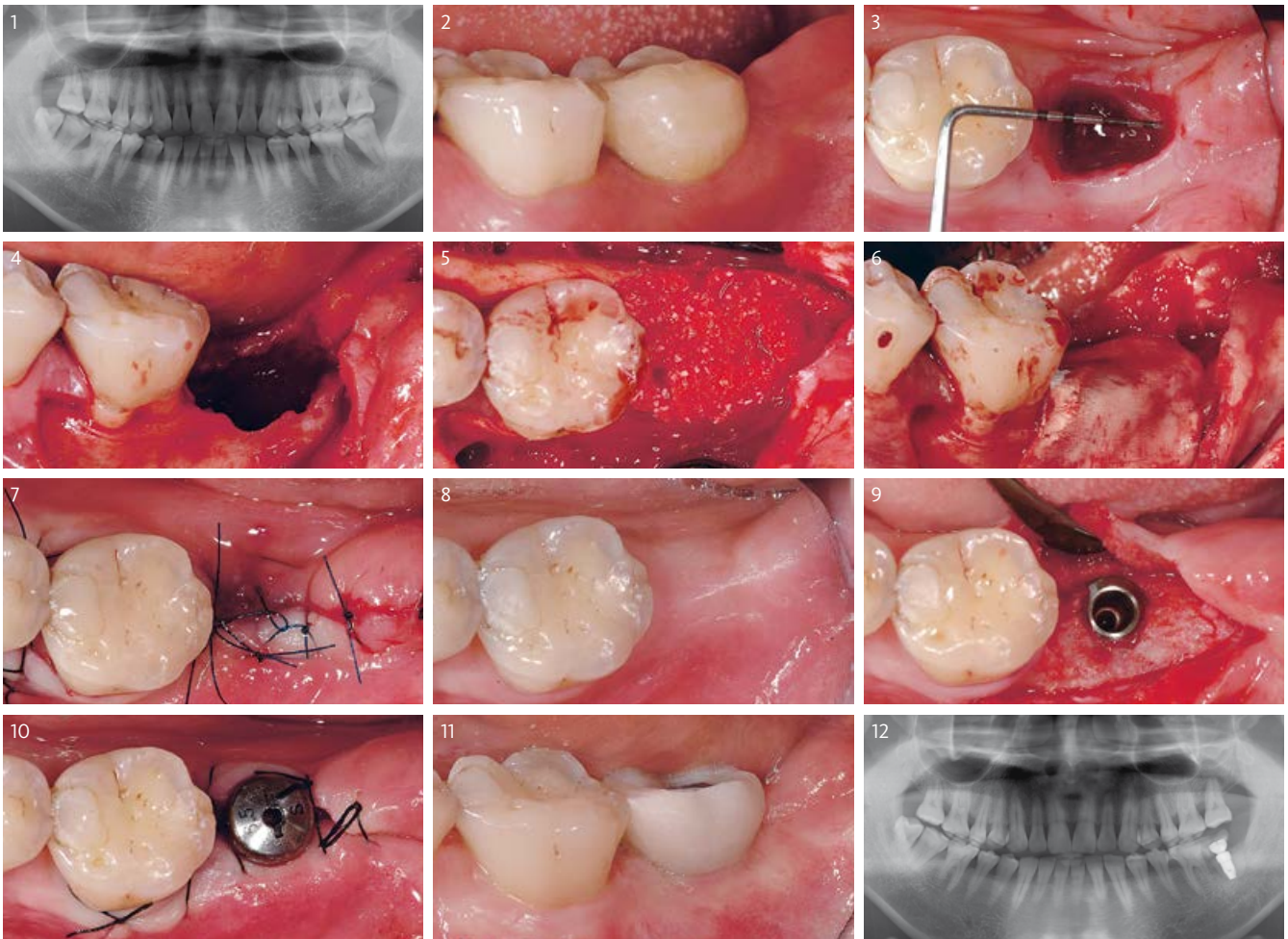

- 1 Radiological status prior to extraction. Initial Smile.
- 2 Starting situation.
- 3 Status following atraumatic extraction of tooth 17.
- 4 A flap is raised.
- 5 Filling of the extraction socket up to the level of the crestal bone level using Geistlich Bio-Oss® Collagen.
- 6 Insertion of the Geistlich Bio-Gide® membrane over the defect .

- 7 Closure of the extraction socket with a mattress suture. Open healing.
- 8 Situation 6 months post-op.
- 9 Newly formed hard tissue. Geistlich Bio-Oss® Collagen is not obvious.
- 10 One stage protocol with healing abutment.
- 11 Provisional prosthesis.
- 12 Radiological view after implantation.

Material selection



Geistlich Combi-Kit Collagen:  
 Geistlich Bio-Oss® Collagen  
 (100 mg)  
 Geistlich Bio-Gide® (16 × 22 mm)

# RIDGE PRESERVATION FOR IMPLANT SUPPORTED BRIDGE

## The patient's risk profile

Aesthetic risk factors	Low risk	Medium risk	High risk
Patient's health	Intact immune system (non-smoker)	Light smoker	Impaired immune system (heavy smoker)
Patient's aesthetic requirements	Low	Medium	High
Height of the smile line	Low	Medium	High
Gingival biotype	Thick "low scalloped"	Medium "medium scalloped"	Thin "high scalloped"
Shape of dental crowns	Rectangular		Triangular
Infections at implantation site	None	Chronic	Acute
Bone height at adjacent tooth	≤ 5 mm from contact point	5.5–6.5 mm from contact point	≥ 7 mm from contact point
Restorative status of adjacent tooth	Intact		Restored
Width of tooth gap	1 tooth (≥ 7 mm)	1 tooth (< 7mm)	2 teeth or more
Soft-tissue anatomy	Intact		Defective
Bone anatomy of the alveolar ridge	No defect	Horizontal defect	Vertical defect

## Quintessence

### Objectives

- › Maintain alveolar contour, which is a combination of hard and soft tissue under pontics.

### Conclusions

- › Geistlich Mucograft® prevents particulate graft from leaking out of the socket before being incorporated into healed tissue.
- › Alveolar contour was largely maintained with Geistlich Mucograft® and Geistlich Bio-Oss®.





“This treatment is ideal for extraction sockets to preserve aesthetic contours when there are limited bony defects.”

Case documentation



- 1 Radiographic findings prior to implant placement in teeth 12 and 22.
- 2 Clinical initial situation prior to implant placement in teeth 12 and 22.
- 3 Maxillary central incisors scheduled for extraction due to recurrent endodontic infections 2 months after implant placement in lateral incisors.
- 4 Extraction sockets grafted with Geistlich Bio-Oss®. The bone substitute fills the socket up to slightly above the bone crest.
- 5 Geistlich Mucograft® is placed over the occlusal surfaces as a socket seal.
- 6 Provisional restoration.
- 7 Provisional restoration contoured to maintain Geistlich Mucograft® in place, taking care not to compress the grafted site.
- 8 Vascularisation and integration of Geistlich Mucograft® after two weeks.
- 9 Clinical situation 1 month post-op.
- 10 Occlusal view at 9 months with the final restoration (11 months after teeth extraction).
- 11 Buccal view at 9 months with the final restoration (11 months after teeth extraction).
- 12 Radiograph showing integration of the graft material in the sockets. Final restoration in place.

Material selection



Geistlich Bio-Oss® small granules (0.25–1 mm)  
 Geistlich Mucograft® (15 × 20 mm punch 8 mm diameter)

# RIDGE PRESERVATION IN MULTIPLE EXTRACTION SOCKETS

## The patient's risk profile

Aesthetic risk factors	Low risk	Medium risk	High risk
Patient's health	Intact immune system (non-smoker)	Light smoker	Impaired immune system (heavy smoker)
Patient's aesthetic requirements	Low	Medium	High
Height of the smile line	Low	Medium	High
Gingival biotype	Thick "low scalloped"	Medium "medium scalloped"	Thin "high scalloped"
Shape of dental crowns	Rectangular		Triangular
Infections at implantation site	None	Chronic	Acute
Bone height at adjacent tooth	≤ 5 mm from contact point	5.5–6.5 mm from contact point	≥ 7 mm from contact point
Restorative status of adjacent tooth	Intact		Restored
Width of tooth gap	1 tooth (≥ 7 mm)	1 tooth (< 7mm)	2 teeth or more
Soft-tissue anatomy	Intact		Defective
Bone anatomy of the alveolar ridge	No defect	Horizontal defect	Vertical defect

## Quintessence

### Objectives

- › Ridge profile maintenance under full arch bridge.
- › Flapless procedure.

### Conclusions


- › Good and quick soft-tissue healing during the early healing phase.
- › Bone volume has been largely preserved with a minimally invasive approach.





“In complex cases, I don’t want to experiment with materials. So I took here the proven Geistlich Biomaterials.”

Case documentation

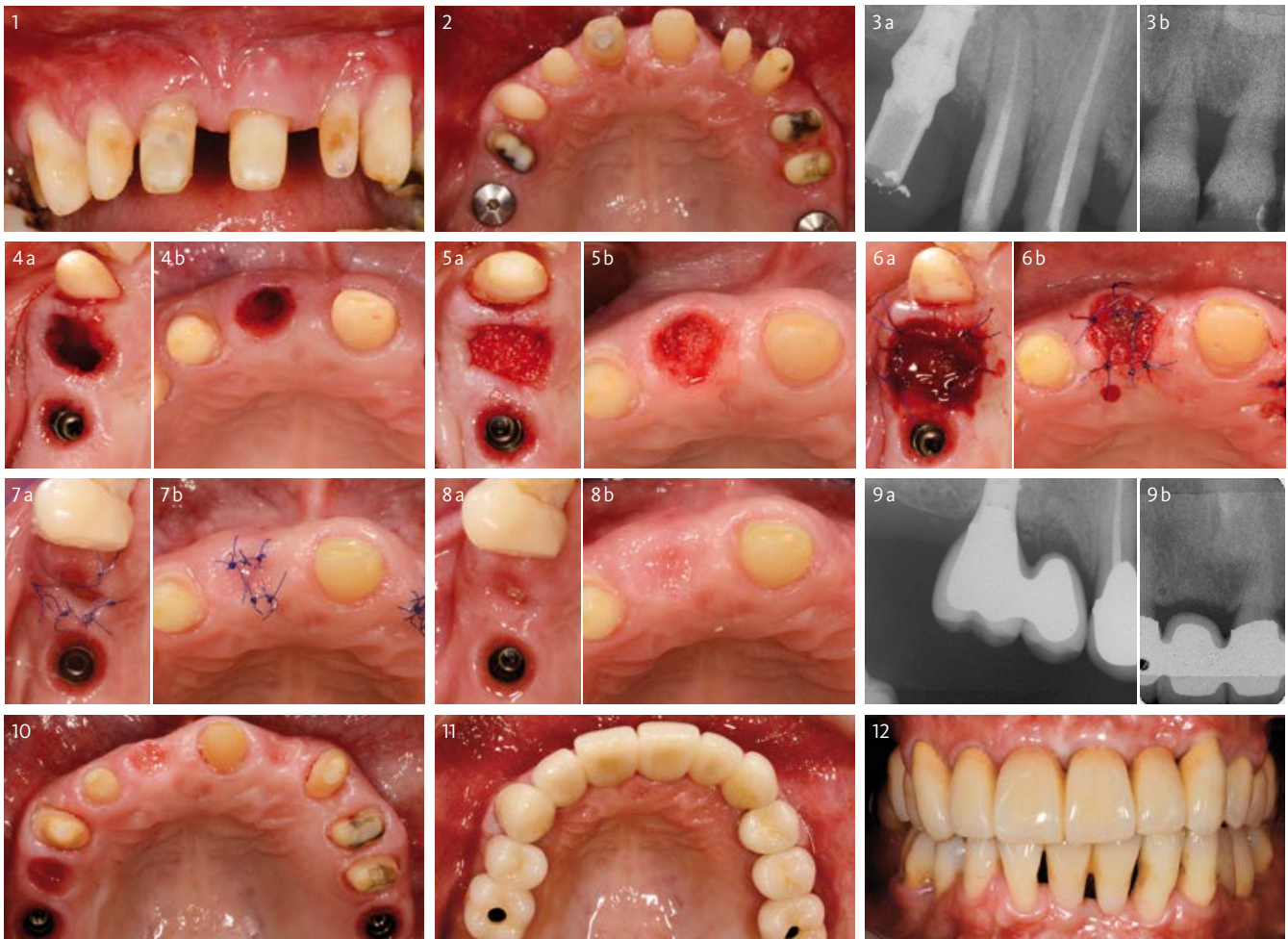

- 1 Initial situation before extraction of teeth 11 and 14.
- 2 Occlusal clinical view showing the ridge profile.
- 3 X-ray findings prior to extraction of teeth a) 14 and b) 11.
- 4 Empty extraction sockets of teeth a) 14 and b) 11.
- 5 Extraction sockets filled with Geistlich Bio-Oss® Collagen.
- 6 Geistlich Mucograft® Seal adapts well to the defects and is sutured with single interrupted sutures.
- 7 Occlusal view before removal of sutures, 1 week after teeth extraction.

- 8 Occlusal view shows nice early healing of the soft-tissues, 1 week post-extraction.
- 9 X-ray findings 12 months post-extraction. Region a) 14 and b) 11.
- 10 Clinical situation of the conditioned soft tissues 12 months post-extraction.
- 11 Final restoration 12 months after extraction (occlusal).
- 12 Final restoration 12 months after extraction (buccal).

Material selection



Geistlich Bio-Oss® Collagen  
(100 mg)  
Geistlich Mucograft® Seal  
(8 mm diameter)

# TECHNICAL GUIDELINES



## Geistlich Bio-Oss® Collagen

- › can be applied both dry, as well as moistened with saline solution or blood.
- › can be cut to size and carefully introduced into the socket with a forceps.
- › can be packed into the socket with a bone graft plugger (or similar), taking care not to compress it too strongly.



## Geistlich Bio-Gide®

- › should be cut dry.
- › should be applied dry with the “UP” mark facing the oral cavity.
- › can be applied inside the alveolus on the defect area or alternatively be inserted between the periosteum and the soft tissue.
- › should be folded over the filled alveolus and adapted under the sulcus.
- › can be fixed by suturing the de-epithelialised soft tissue over the exposed membrane, (e.g. with single sutures) or completely submerged under a flap.



## Geistlich Mucograft® Seal<sup>1</sup>

- › has to be used with an alveolar filling material (e.g. Geistlich Bio-Oss® Collagen).
- › should be applied after de-epithelialisation of the adjoining soft-tissue margins.
- › should be adapted to the defect size and applied dry.
- › has to be applied with the spongy framework (marked with grooves) facing towards the extraction socket.
- › should be sutured with non-resorbable suture and not glued.
- › should be sutured with single-interrupted sutures (recommended: 5.0 or 6.0), double interrupted sutures or cross sutures (recommended: 5.0), depending on the defect.
- › should be tension-free and closely adapted to the de-epithelialised marginal soft-tissue border.



## Reference

- <sup>1</sup> Adapted from Geistlich Mucograft® Seal Advisory Board Meeting Report 2013. Data on file, Geistlich Pharma AG, Wolhusen, Switzerland.



# PRODUCT RANGE



## Geistlich Bio-Oss®

Small granules (0.25–1 mm) | Quantities: 0.25 g, 0.5 g, 2.0 g (1 g ≈ 2.05 cm<sup>3</sup>)  
Large granules (1–2 mm) | Quantities: 0.5 g, 2.0 g (1 g ≈ 3.13 cm<sup>3</sup>)

The small Geistlich Bio-Oss® granules are recommended for smaller 1–2 socket defects and for contouring autogenous block grafts. The large Geistlich Bio-Oss® granules enable improved regeneration over large distances and provide enough space for the ingrowing bone.



## Geistlich Bio-Oss Pen®

Small granules (0.25–1 mm) | Quantities: 0.25 g ≈ 0.5 cc, 0.5 g ≈ 1.0 cc  
Large granules (1–2 mm) | Quantity: 0.5 g ≈ 1.5 cc


Geistlich Bio-Oss® granules are available in an applicator. It allows the bone substitute material to be applied faster and more precisely to the surgical site. Geistlich Bio-Oss Pen® is available containing both the small granules and the large granules.



## Geistlich Bio-Oss® Collagen

Geistlich Bio-Oss® (small granules) + 10% collagen (porcine)  
Sizes: 100 mg (0.2–0.3 cm<sup>3</sup>), 250 mg (0.4–0.5 cm<sup>3</sup>), 500 mg (0.9–1.1 cm<sup>3</sup>)

Geistlich Bio-Oss® Collagen is indicated for use in periodontal defects and extraction sockets. Through the addition of collagen, Geistlich Bio-Oss® Collagen can be tailored to the morphology of the defect and is particularly easy to apply.

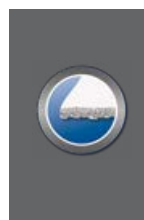

## Geistlich Bio-Gide®

Resorbable bilayer membrane  
Sizes: 25 × 25 mm, 30 × 40 mm

Geistlich Bio-Gide® consists of porcine collagen and has a bilayer structure – a rough side that faces the bone and a smooth side that faces the soft tissue. Geistlich Bio-Gide® is easy to handle: it can be positioned easily, adheres well to the defect, and is resistant to tension and tearing.



## Geistlich Combi-Kit Collagen

Geistlich Bio-Oss® Collagen 100 mg  
+ Geistlich Bio-Gide® 16 × 22 mm


When used in combination, the system has optimised properties for Ridge Preservation and minor bone augmentations according to the GBR principle.



## Geistlich Mucograft® Seal

Collagen matrix  
Size: 8 mm diameter

Geistlich Mucograft® Seal consists of a compact structure that gives stability while allowing open healing, and a spongy structure that supports blood clot stabilisation and ingrowth of soft-tissue cells.

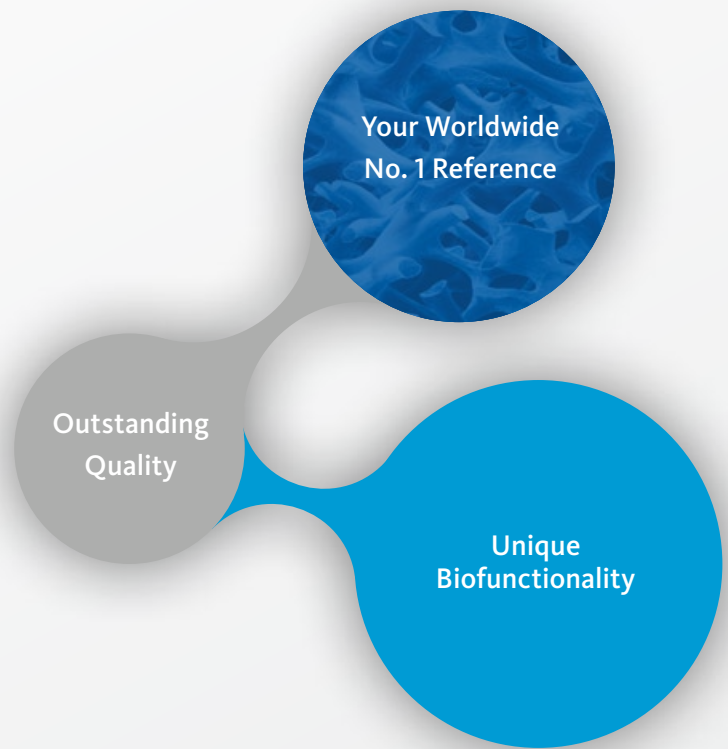

## Your Worldwide No. 1 Reference

Geistlich Biomaterials is constantly working to offer you solutions for easy, predictable and successful management and regeneration of extraction sockets. The company's own research departments along with global experts develop the product portfolio, and try new techniques and applications for existing products. In more than 15 worldwide Round Table Meetings\*, expert clinicians and Geistlich Biomaterials cooperate on the aim of promoting discussion and evolving a consensus on the treatment concepts for extraction sockets. These Round Table Meetings also help to define what is the current published clinical evidence and where research still needs to be done.

## Outstanding Quality

Quality and safety are high priorities at Geistlich Pharma. At Geistlich Pharma everything is done under one roof: from the selection and control of the raw material to production and storage until dispatch, and all steps are taken seamlessly and meet the company's high standards of quality and safety.

## Unique Biofunctionality

The excellent results of Ridge Preservation with Geistlich Biomaterials are largely due to their unsurpassed biofunctionality: Geistlich Bio-Oss® with its porous structure<sup>1</sup> serves as guide rail for the in-growing blood vessels<sup>2</sup> and integrates into newly formed bone<sup>3</sup>, while the unique bilayer Geistlich Bio-Gide® shields the young bone from the surrounding connective tissue cells and supports wound healing<sup>4</sup> and early vascularisation<sup>5</sup>. The 3-dimensional matrix of Geistlich Mucograft® Seal facilitates soft-tissue cells ingrowth<sup>6</sup> and may enhance early wound healing<sup>7</sup>.

Clinically relevant:

- › Geistlich Biomaterials are perfectly suited to combined use for treatment of extraction sockets.
- › Geistlich Bio-Oss® Collagen combined with Geistlich Bio-Gide® preserves up to 93% of the ridge width<sup>8,9</sup> and they promote more new bone formation vs. no membrane<sup>10</sup>
- › Geistlich Bio-Oss® Collagen combined with Geistlich Mucograft® Seal increases preserved bone volume when compared to spontaneous healing<sup>11</sup>

### References

- 1 Weibrich G et al., Mund Kiefer Gesichtschirurg 4, 2000; 148–152.
- 2 Degidi M et al., Oral Dis. 2006 Sep; 12(5): 469–475.
- 3 Artzi Z, et al. J Periodontol. 2001 Feb;72(2):152-9.
- 4 Becker J et al., Clin. Oral Implants Res. 2009; 20(7): 742–93.
- 5 Rothamel D et al., Clin. Oral Implants Res. 2005;16:369–378.
- 6 Ghanaati S, et al. Biomed Mater. 2011 Feb;6(1):015010.
- 7 Thoma DS, et al. J Clin Periodontol. 2012 Feb;39(2):157-65.
- 8 Cardaropoli D, et al. Int J Periodontics Restorative Dent. 2012 Aug;32(4):421-30.
- 9 Cardaropoli D, et al. Int J Periodontics Restorative Dent. 2014 Mar-Apr;34(2):211-7.
- 10 Perelman-Karmon et al. Int J Periodontics Restorative Dent. 2012 Aug;32(4):459-65.
- 11 Jung RE, et al. J Clin Periodontol. 2013 Jan;40(1):90-8.

\* Data on file (Wolhusen, Switzerland): Austria, Baltics, Belgium, Brazil, France, Germany, Greece, Holland, Korea, Nordics, Poland, Russia, Spain (2009), Spain/Portugal (2014, 2015), Switzerland (2009, 2011, 2013), UK.

GET THE NO. 1  
INTO YOUR PRACTICE



**Manufacturer**

© Geistlich Pharma AG  
Business Unit Biomaterials  
Bahnhofstrasse 40  
CH-6110 Wolhusen  
Phone +41 41 492 56 30  
Fax +41 41 492 56 39  
[www.geistlich-biomaterials.com](http://www.geistlich-biomaterials.com)

**Subsidiary Great Britain,  
Ireland**

Geistlich Biomaterials  
Geistlich Sons Limited  
1<sup>st</sup> Floor, Thorley House  
Bailey Lane  
Manchester Airport  
GB-Manchester M90 4AB  
Phone +44 161 490 2038  
Fax +44 161 498 6988  
[www.geistlich.co.uk](http://www.geistlich.co.uk)

More details about our  
distribution partners:  
[www.geistlich-biomaterials.com](http://www.geistlich-biomaterials.com)